

Shaping high-end imaging for your research



ZEISS LSM 980 with Airyscan 2

Your next generation confocal for fast and gentle multiplex imaging

Your Next Generation Confocal for Fast and Gentle Multiplex Imaging

› In Brief

› The Advantages

› The Applications

› The System

› Technology and Details


› Service

Life sciences research can be demanding, and if you are involved in neuroscience, cancer research or other cell- or organism-based disciplines, you'll often need microscopy data for your work. Emerging technologies such as CRISPR/Cas open up innovative ways of thinking and allow you to ask altogether new scientific questions, deeply affecting your imaging experiments. To monitor life as undisturbed as possible requires low labeling density for your biological models—for example, 3D cell culture, spheroids, organoids or even whole organisms—and this calls for 3D imaging that combines optical sectioning with low phototoxicity and high speed. Then there are the repeated experiment runs it takes to get statistically-valid data for your conclusions: it soon becomes apparent you will also need high throughput.

LSM 980 with Airyscan 2 is the ideal platform for your confocal 4D imaging, optimized for simultaneous spectral detection of multiple weak labels up to 900 nm emission with the highest light efficiency.

Add Airyscan 2 with its Multiplex mode to get more imaging options for your experiments. You can now choose the perfect setup to gently image larger fields of view with super-resolution in shorter acquisition times than ever before.

A number of software helpers will optimize your workflow and support efficient acquisition and data management. With ZEN Connect you can document and share all details of your experiments. You'll always keep the context as you combine overview images, ROIs and additional data, even across imaging modalities.

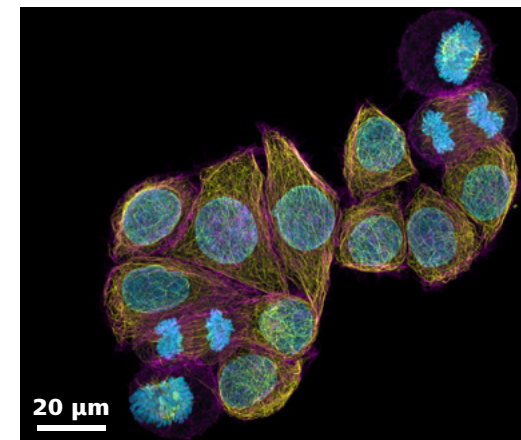

HeLa cells stained for DNA (blue, Hoechst 44432), microtubules (yellow, anti-tubulin Alexa 488) and F-actin (magenta, phalloidin Abberior STAR Red). Imaged with ZEISS Airyscan 2 in Multiplex mode for efficient super-resolution imaging of a large field of view. Courtesy of A. Politi, J. Jakobi and P. Lenart, MPI for Biophysical Chemistry, Göttingen, Germany.

See for yourself how the Multiplex mode for Airyscan 2 gives you better data faster than ever before. Book a hands-on demonstration in one of our ZEISS Microscopy Labs now.


>> www.zeiss.com/lsm980

Simpler. More Intelligent. More Integrated.

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

Get Better Data Faster


Use the Multiplex mode for Airyscan 2 and get more information in less time. Smart illumination and detection schemes let you image your most challenging three-dimensional samples with high framerate rates beyond the diffraction limit and still treat your sensitive samples gently. By combining the full flexibility of a point scanning confocal with the speed and gentleness of the sensitive Airyscan area detector, it's now possible to answer your scientific questions up to ten times faster with super-resolution.

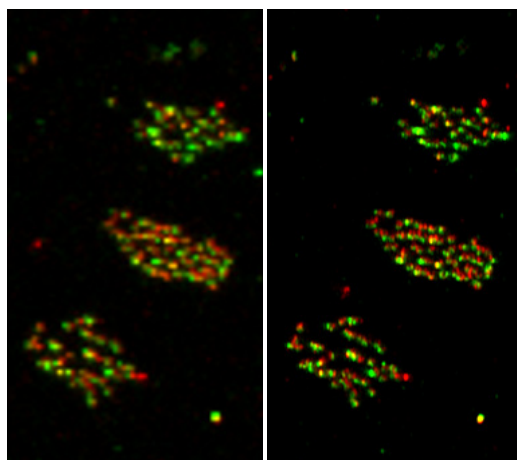

Basal bodies (red) and basal feet (green) of ependymal cilia labelled by immunohistochemistry. Left side: imaged in conventional confocal mode; right side: imaging with the same frame time using the Multiplex mode for Airyscan 2 clearly reveals the orientation of the cilia. Courtesy of S. Kapoor, Max Planck Institute for Biophysical Chemistry, Göttingen, Germany.

Image with More Sensitivity


LSM 980 brings you the best of two worlds to image your most challenging samples. You get the light efficient beam path of the LSM 9 family with up to 36 simultaneous channels for full spectral flexibility up to the near infrared (NIR) range. Plus, when you combine it with Airyscan 2, this revolutionary area detector extracts even more information from your sample in less time. This lets you image faint signals with 4 – 8 times higher signal-to-noise ratio (SNR). You don't need to close a pinhole to get super-resolution, which makes your 3D imaging even more light efficient.

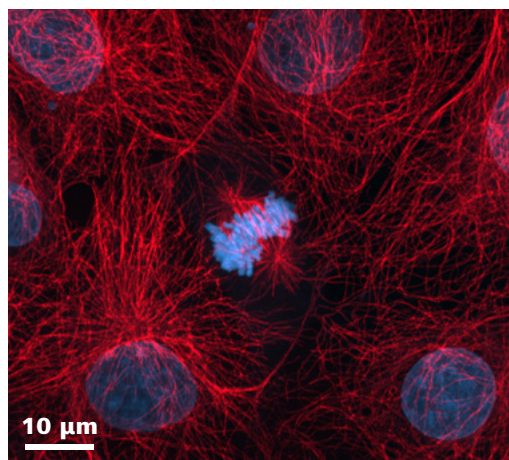

Cos-7 cells, DAPI (blue) and Anti-Tubulin Alexa 700 (red). Alexa 700 imaged with NIR detector. Sample courtesy of Urs Ziegler and Jana Doehner, University of Zurich, ZMB, Switzerland.

Increase Your Productivity


It's never been easier to set up complex confocal live cell imaging experiments. ZEN microscopy software puts a wealth of helpers at your command to achieve reproducible results in the shortest possible time. AI Sample Finder helps you quickly find regions of interest, leaving more time for experiments. Smart Setup supports you in applying best imaging settings for your fluorescent labels. Direct Processing enables parallel acquisition and data processing. ZEN Connect keeps you on top of everything, both during imaging and later when sharing the whole story of your experiment. It's easy to overlay and organize images from any source.

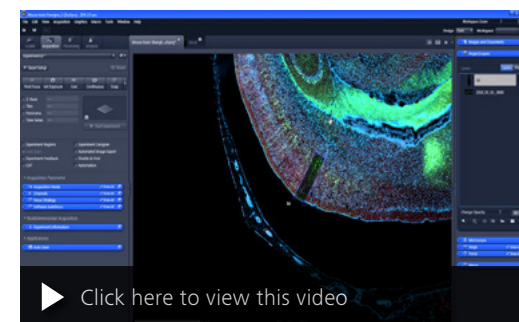

See how ZEN Connect helps to always keep your context while imaging. From acquiring an overview image, to defining ROI's, and even when changing between different imaging systems. You save time and always stay on top of things.

Your Insight into the Technology Behind It

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

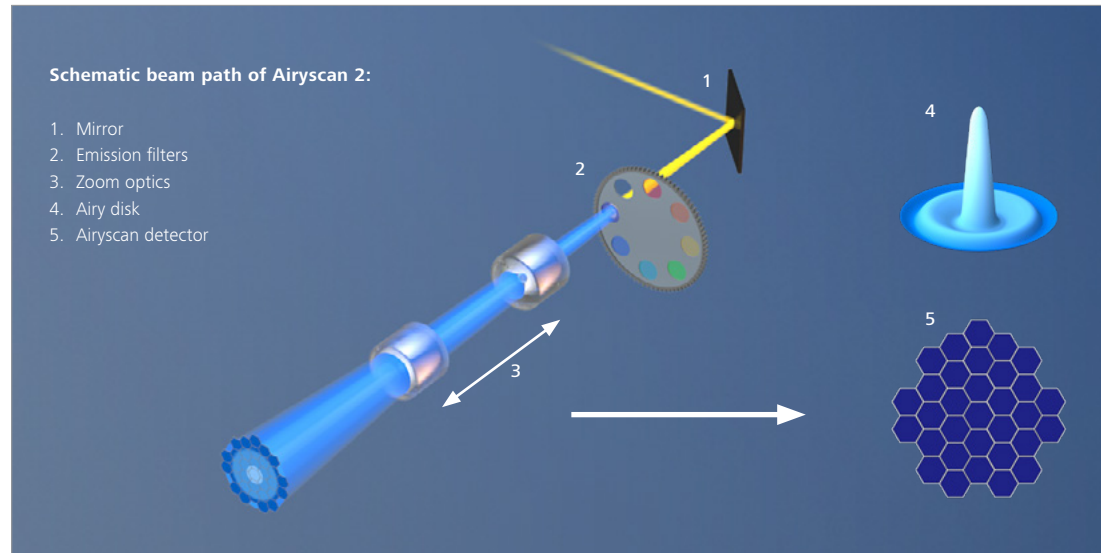



The Airyscan Principle

Classic confocal laser scanning microscopes use point illumination to scan the sample sequentially. The microscope optics transform each point to an extended Airy disk (Airy pattern). A pinhole then spatially limits this Airy disk to block out-of-focus light from reaching the detector. Closing the pinhole gives higher resolution, but at the price of detecting fewer photons – and these photons cannot be brought back by e.g. deconvolution.

Airyscan 2 is an area detector with 32 concentrically arranged detection elements. This allows you to acquire more of the Airy disk at once. The confocal pinhole itself remains open and does not block light, thus more photons are collected. This produces much greater light efficiency while imaging. Airyscan 2 gives you a unique combination of gentle super-resolution imaging and high sensitivity.

For further information on the Airyscan principle please refer to:

<https://zeiss.ly/airyscan-principle>

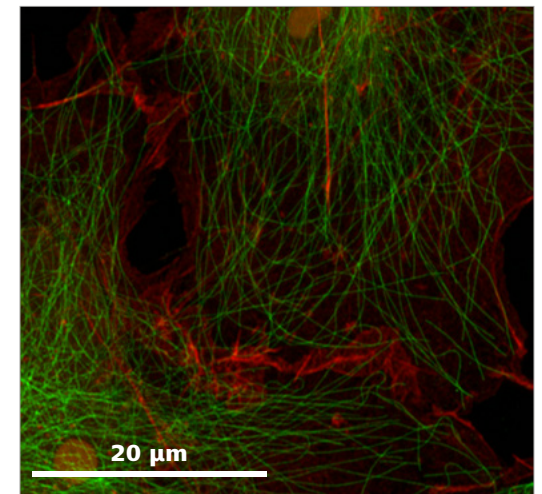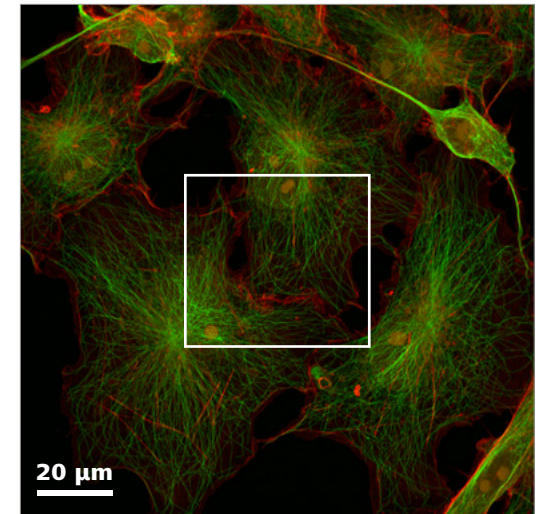

Comparing the field of view you can image at super-resolution in the same time using Airyscan SR (bottom) and Multiplex mode (top). Cos-7 cells with labelled microtubules (alpha-tubulin 488, green) and actin (phalloidin 647, red).

Your Insight into the Technology Behind It


- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

The Multiplex Mode for Airyscan 2


Do you want to image large fields of view and whole sample volumes in the shortest time possible? And do you want to image with superb image quality at the same time? The LSM 9 family with Airyscan 2 from ZEISS now gives you more options to fit imaging speeds and resolution to your experimental needs. You combine an area detector with smart illumination and readout schemes, which let you choose from different parallelization options.

The Multiplex mode uses knowledge about the shape of the excitation laser spot and the location of single area detector elements to extract more spatial information, even during parallel pixel readout. This allows to take bigger steps when sweeping the excitation laser over the field of view, improving your achievable acquisition speeds. In fact, the high amount of spatial information captured in the pinhole plane allows to reconstruct a final image with better resolution than the acquisition sampling. Airyscan 2 in Multiplex mode can acquire up to four super-resolution image lines with high SNR in a single sweep. Your LSM 980 with Airyscan 2 allows to stretch the excitation laser spot to image eight lines in parallel. Use this speed advantage for ultrafast time series of single slices, for rapid tiling of large areas or for fast volumetric time-lapse imaging.

Choose the Perfect Airyscan Mode for Your Experiment

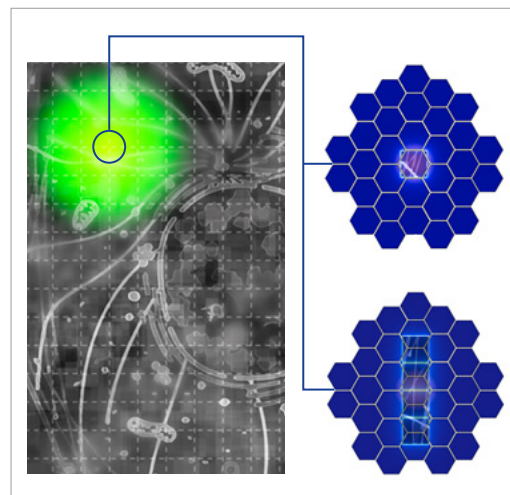

For each illumination position, Airyscan SR mode generates one super-resolution image pixel. The spatial information provided by Airyscan 2 in Multiplex SR-4Y allows to scan 4 super-resolution image lines in a single sweep.

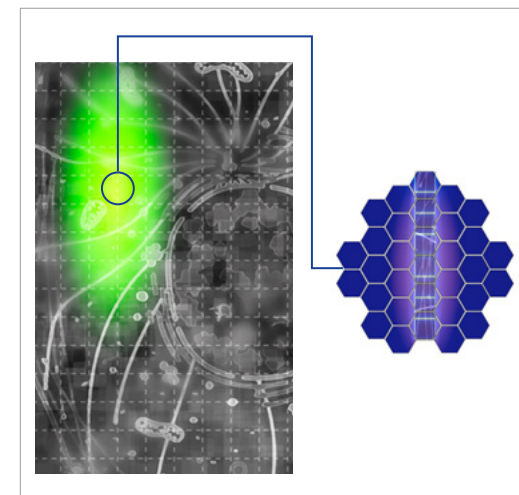

For Airyscan Multiplex SR-8Y and CO-8Y the illumination laser spot is vertically elongated which allows to capture 8 image pixels for each illumination position. Sampling can be done in super-resolution (SR) or confocal (CO) resolution, depending on your experiment.

LSM 980 with Airyscan 2


	Airyscan SR	Multiplex SR-4Y	Multiplex SR-8Y	Multiplex CO-8Y
Parallelization	1	4	8	8
Resolution	120/120	140/140	120/160	180/220
FPS at max FOV	0.2 (Zoom 1.7)	1.0 (Zoom 1)	2.0 (Zoom 1)	9.6 (Zoom 1)
FPS at 512 × 512 pixels	4.7	25	47.5	34.4
Antibody labeling, fine structures	+++++	++++	+++	++
Antibody labeling, tiling	++	++++	+++++	+++
Live cell imaging	++	+++	++++	+++++

Your Insight into the Technology Behind It

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

Beam path of LSM 980 with Airyscan 2

Emission light travels through the Twin Gate main dichroic beam splitter with its very efficient laser suppression to deliver supreme contrast. Then, at the secondary beam splitter, the emission light either travels to the internal spectral detection unit (Quasar) with 3, 6, or 34 channels and to the NIR detection unit. Or, light is sent to the Airyscan 2 area detector with GaAsP technology.

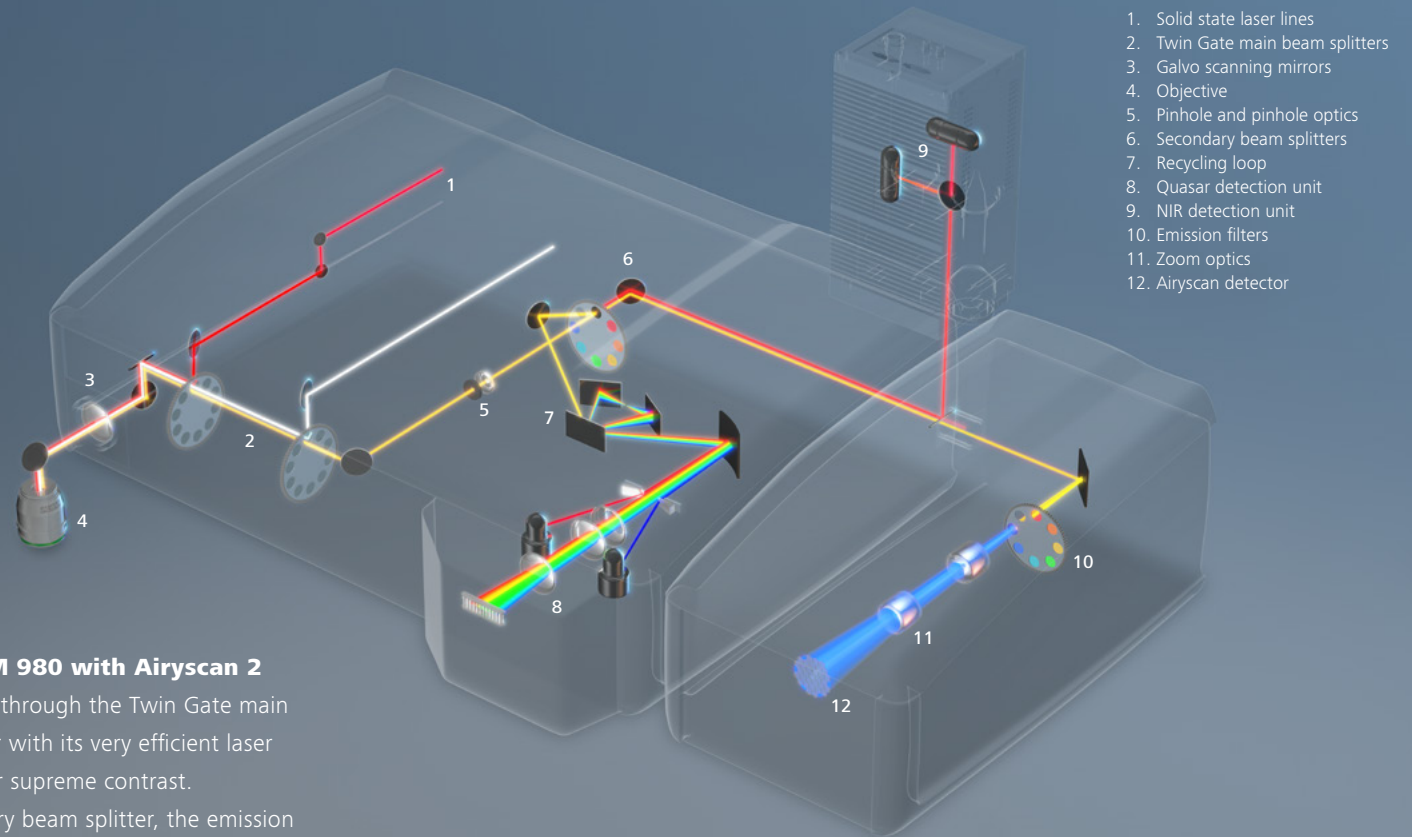

▶ [Click here to view this video](#)


Your Insight into the Technology Behind It

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service



LSM 980 brings a great deal of freedom to your experimental setup. The beam path design and every single component are optimized to deliver the highest sensitivity and spectral flexibility for your experiments.

Sensitivity

Sensitivity is key to visualize the lowest signal and resolve all structures. The LSM 980 beam path design ensures that you don't need to worry about this part of your experiments. The linear galvo scanners ensure that each image pixel gets identical time to contribute to the image and more than 80% of each frame time is used to collect your precious signal. The low angle Twin Gate beam splitter directs the excitation laser light towards your sample and keeps it efficiently separated from your emission signal. You can even cover the detection range over laser lines to collect all emission light, without worrying about stray light or laser light reflection. For your multi-label experiments, each emission color is captured with the most suited detector technology along the complete wavelength range from 380 nm up to 900 nm. Homogeneous spectral separation is achieved with the holographic grating. The unique recycling loop directs all your photons towards the 3, 6, or 34 channel Quasar detection unit that allows you to define the emission detection bands in nanometer precision with ease. You can use the lowest laser powers even at extremely low expression levels and capture your information by switching GaAsP or GaAs detectors to photon counting mode.

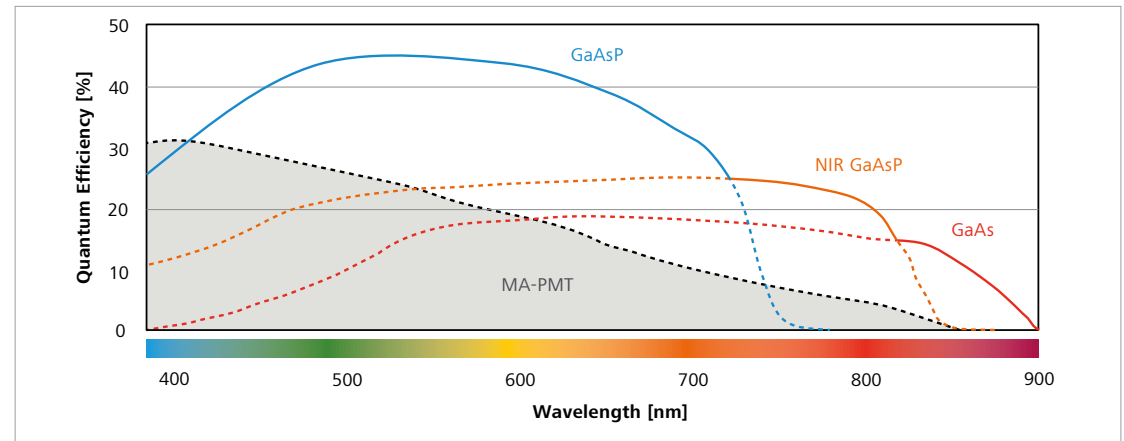

Typical spectral quantum efficiency (QE) of ZEISS LSM 980 detectors

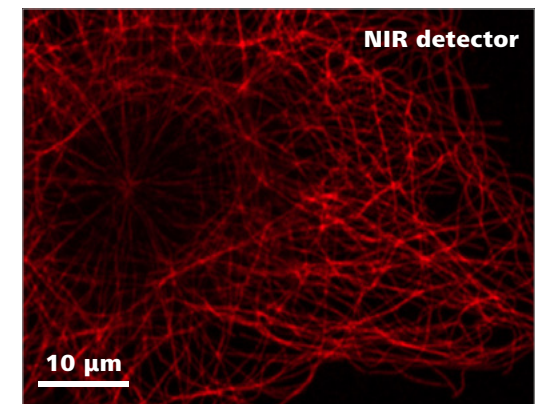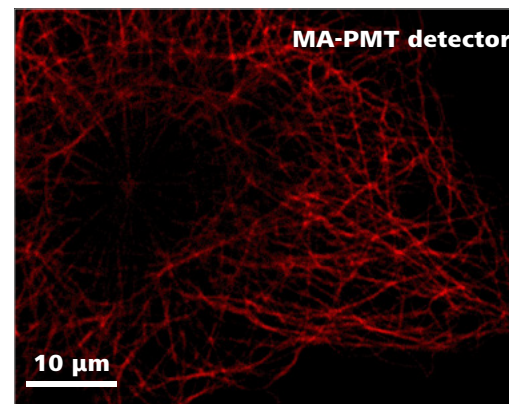


Microtubules of Cos-7 cell (Anti-Tubulin AF700). Comparison of the ZEISS LSM 980 MA-PMT and the ZEISS NIR GaAsP detector; excitation with 639 nm laser at same laser power. Emission range for the MA-PMT is set to 660 – 757 nm, and for the NIR detector is 660 – 900 nm. Sample courtesy of Urs Ziegler and Jana Doehner, University of Zurich, ZMB, Switzerland.

Your Insight into the Technology Behind It

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

Spectral flexibility

Freely select the fluorescent labels from 380 nm to the near infrared (NIR) range that are best for your experiment and combine them as needed. Smart Setup supports you with the best multi-fluorescent image setup by choosing suitable solid-state laser lines for excitation and the detectors with the highest quantum efficiency for each wavelength. Expanding your spectral range into the NIR allows you to use more labels in parallel. Visualize additional structures with more dyes in your multi-color experiments, with the Quasar and NIR detectors efficiently supporting spectral multiplexing experiments. NIR fluorescent labels are less phototoxic for living samples due to the longer wavelength. This allows you to investigate living samples for longer periods of time while limiting the influence of light. Additionally, light of longer wavelength ranges is less scattered by the sample tissue, increasing penetration depth. For any of the advantages you pursue with NIR labels, the dual-channel NIR detector combines two different detector technologies (extended red GaAsP and GaAs) for optimal sensitivity up to 900 nm. To separate even highly overlapping signal or to remove autofluorescence, you can take a Lambda Scan using the complete detection range with up to 36 detectors, keeping illumination and time needed at a minimum.

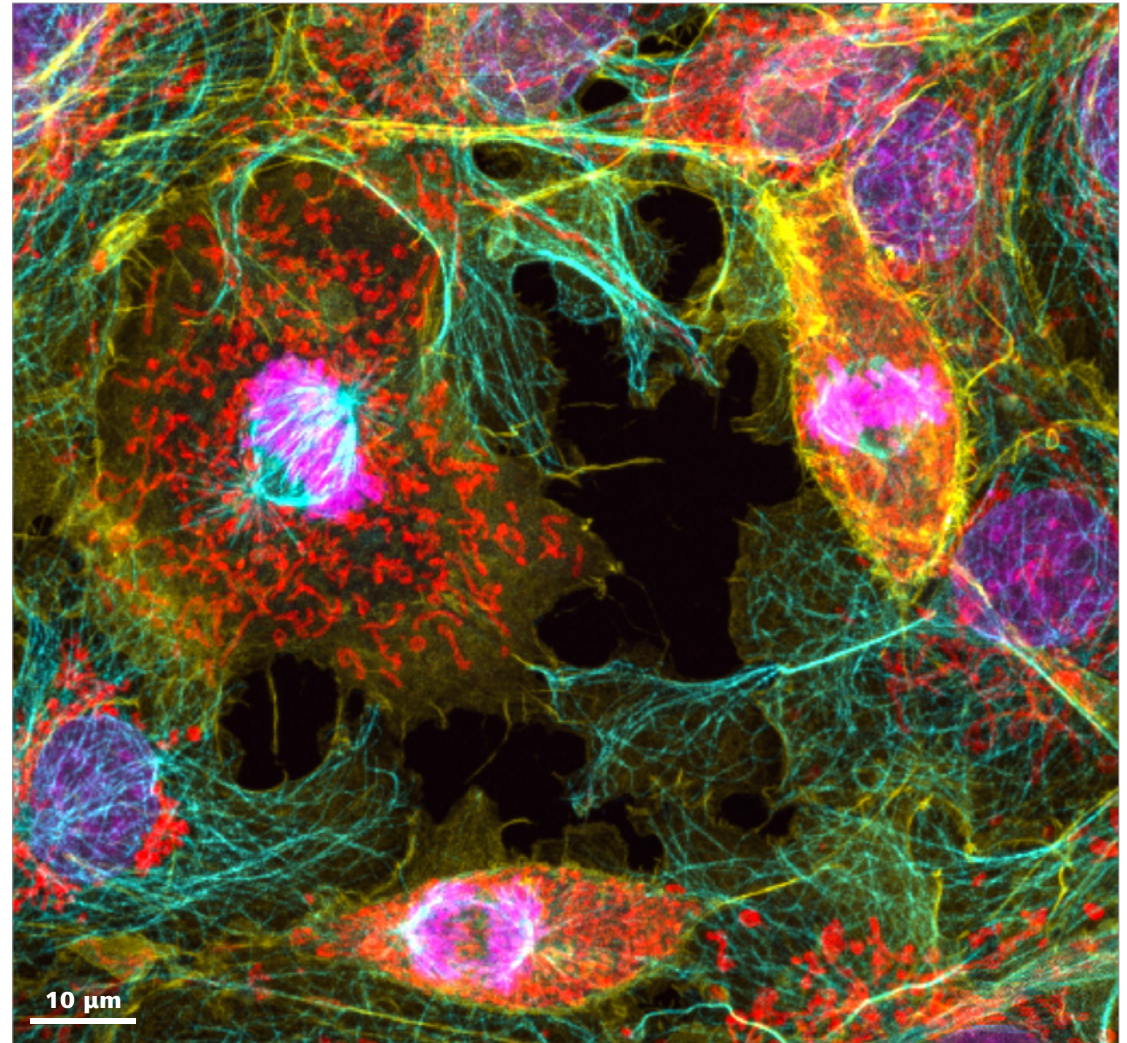

Cos-7 cells, DAPI (magenta), Anti-tubulin Alexa 568 (blue), Actin Phalloidin-OG488 (yellow) and Tom20-Alexa 750 (red). Imaged in Lambda mode across the visible and NIR spectrum. Individual signals separated by Linear Unmixing. Maximum intensity projection of a z-stack. Sample courtesy of Urs Ziegler and Jana Doehner, University of Zurich, ZMB, Switzerland.

Your Insight into the Technology Behind It


- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

Data beyond imaging

The combination of laser point illumination, linear scanning and detectors that are able to capture the signal in photon counting mode makes the LSM 980 more than an imaging device:


- **Raster Image Correlation Spectroscopy (RICS)** can generate a display map of molecule concentration and diffusion coefficients of a complete image frame of e.g., a cell.
- **Fluorescence Correlation Spectroscopy (FCS)** allows a non-invasive insight into molecular concentration and diffusion processes, leading to a deeper understanding of cell functions. To measure on a single molecule basis, you can use single- or multiphoton laser lines and use the full emission range up to 900 nm.
- **Fluorescent Cross Correlation Spectroscopy (FCCS)** allows you to observe molecular interaction by utilizing the 32 channels of the Quasar and perform with up to 7 individual channels.
- **Förster Resonance Energy Transfer (FRET)** is another method to investigate molecular distances or interactions, using sensitized emission or acceptor photobleaching approaches.
- **Fluorescence Recovery after Photobleaching (FRAP)** utilizes any of the laser lines to perform flexible photobleaching experiments. The same principle adheres to photomanipulation experiments in general, for example to investigate intracellular movement. Or follow cellular movement within whole organisms by photoconversion of fluorescent protein labels.
- **Fluorescence Lifetime Imaging Microscopy (FLIM)** uses differences in the fluorescence decay to separate components. It is used for functional imaging utilizing that fluorescence lifetime can be influenced by multiple factors such as ion or oxygen concentration, pH and temperature. FLIM is beneficial for FRET measurements, analyzing proximity and interaction of molecules.*

*available upon request

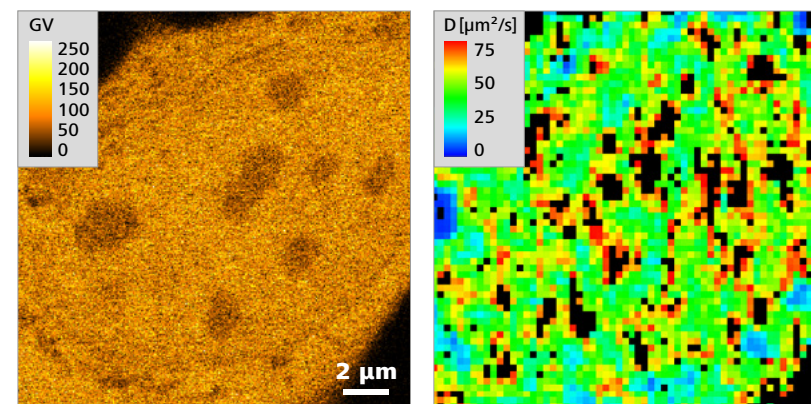

RICS measurement, using U2OS cells expressing monomeric eGFP. The diffusion of the target can be displayed as a map (right) based on the intensity image (left).

Sample courtesy of P. Hemmerich, Leibniz Institute on Aging – Fritz Lipmann Institute, Jena Germany

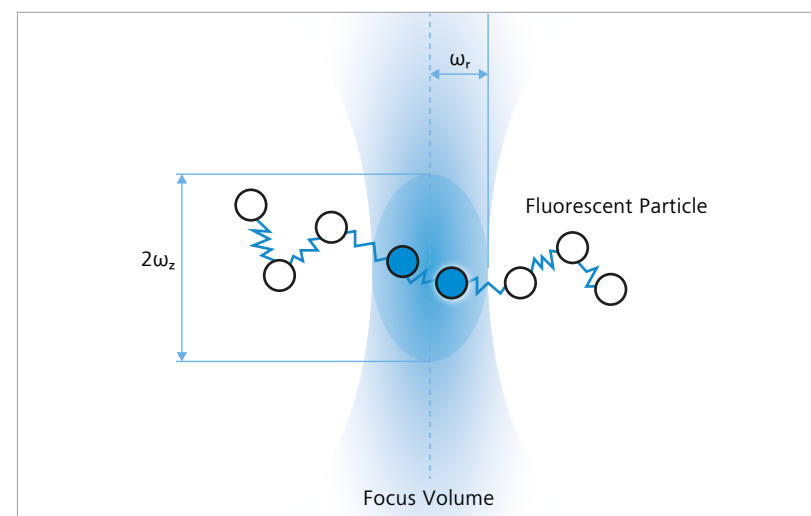

Fluorescence Correlation Spectroscopy (FCS) principle. Trajectory of a fluorescent particle through the detection volume

Your Insight into the Technology Behind It

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

AI Sample Finder: Automated sample identification for efficient imaging

Microscopes are becoming increasingly automated. For sample placement, however, microscope parts such as the condenser arm often have to be moved manually. Focus adjustment and identification of the relevant areas on the sample carrier require additional manual steps.

The AI Sample Finder automates this sequence, eliminating time-consuming manual adjustments and reducing the time to image from minutes to just seconds.


You can access all sample areas directly which allows you starting your experiment faster than ever. The AI Sample Finder greatly improves productivity as you can easily image only those regions containing sample not overlooking potentially important areas.




- After you placed the sample on the loading position, the AI Sample Finder automatically moves it to the objective.



- Without the need of manual sample positioning or focusing, an overview image for fast and convenient navigation is taken within seconds. Composite darkfield illumination creates a high-contrast image even for very low-contrast samples.

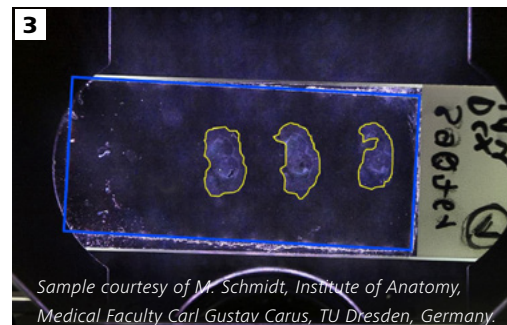

- Intelligent routines automatically identify your sample carrier, regardless if you use a petri dish, a chamber slide, or a multiwell plate. Carrier properties are automatically transferred to the software, eliminating manual settings.

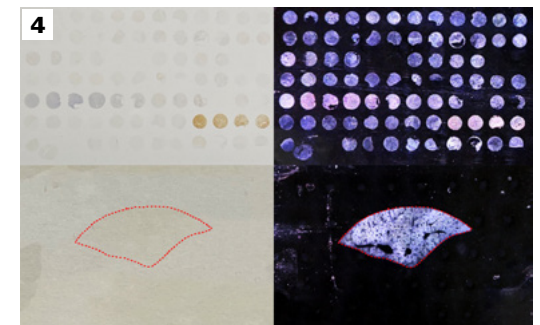

- Your samples are reliably identified. Deep Learning algorithms precisely detect even unusual sample regions. You can navigate and access all sample areas directly which allows you starting your experiment faster than ever.


Expand Your Possibilities

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

Acquire Reproducible Data with Ease

With all its various aspects and workflows, your research leaves you with no time to waste. That's why ZEN microscopy software was created—to make your confocal imaging both efficient and enjoyable. ZEN – ZEISS Efficient Navigation – is the only user interface you will ever see on all imaging systems from ZEISS. This familiar and easy-to-learn interface will help you get reproducible results in the shortest possible time.

Use Smart Setup to select your dyes and ZEN will automatically apply all necessary settings for all LSM imaging modalities. The integrated database with spectral data for more than 500 dyes helps you make an informed decision about your imaging options. You can always save imaging configurations or even whole experiments to reproduce settings quickly. The Reuse function allows you to extract and load imaging settings from the existing images. You will be amazed how easy imaging becomes when the AI Sample Finder automatically detects the sample carrier, adjusts the focus, and finds the sample regions relevant for your experiment. It takes less time to illuminate your sample and leaves you more of the precious time you've booked on the system for imaging. In addition, you can use the overview image to document all steps of your experiment and load it in ZEN Connect to combine with other multi-modal data or aspects of your sample.

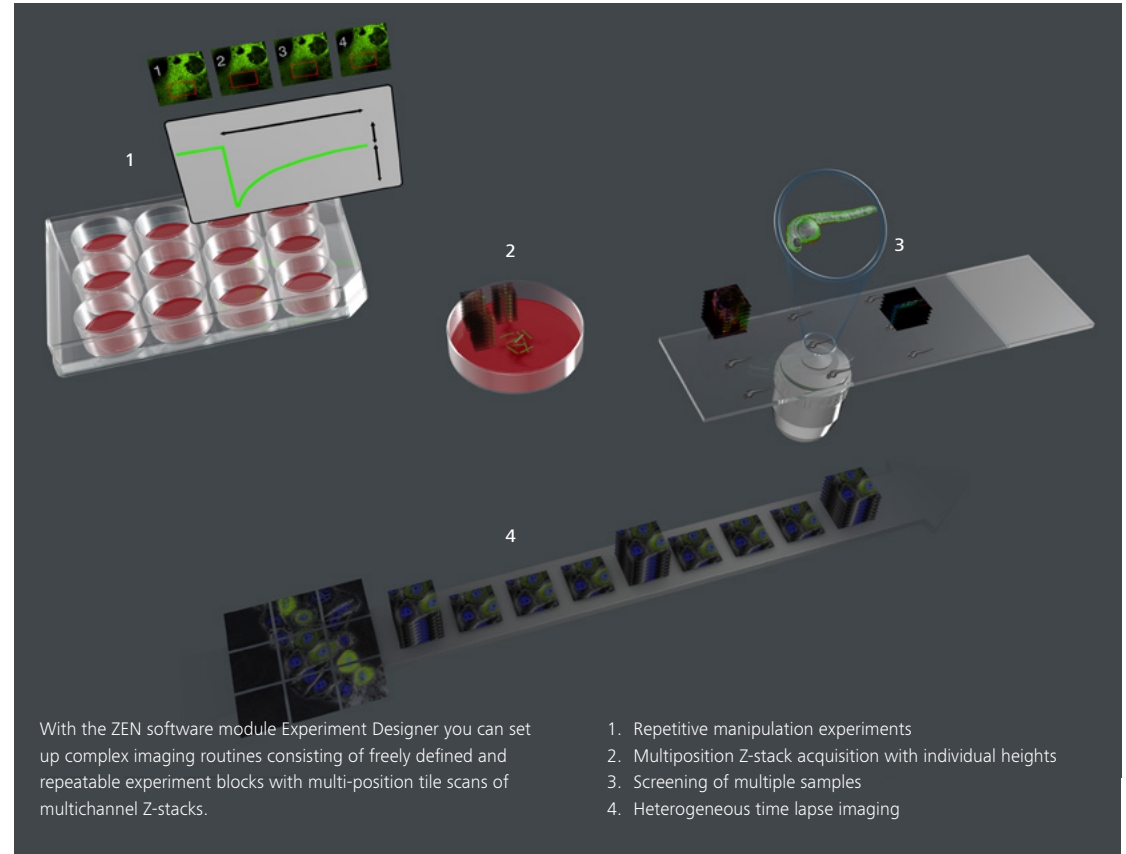

Sometimes your scientific questions will require complex acquisition strategies. Statistical analysis might call for repetitive imaging of a large number of samples with the same or even differing imaging conditions. Experiment Designer is a powerful yet easy-to-use module that images multiple regions with all imaging modalities of your LSM 980.

It gives you access to a number of hardware and software options which will always keep your sample in focus, even during the most demanding long-term time-lapse experiments. You can even view and save your valuable data during acquisition sessions to assess, analyze and react immediately.


Expand Your Possibilities

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

See More Details

Sometimes you need to see and assess your multimodal images during acquisition in order to plan your next steps. ZEN microscopy software gives you multiple options. You can sit at your connected computer to start the Direct Processing function for processing your Airyscan images during acquisition.

However, confocal imaging is only one part of the big picture, and you may need data from additional imaging modalities to complement the view on your sample. ZEN Connect can bring information from all your experiments together. Keep the context of your data by collecting all images of one experiment session in a single project in which you can combine overview and detailed high-resolution images, all perfectly aligned. Once you have created a project, you can always add and align content from any other imaging source, be it ZEISS, non-ZEISS or even cartoons and analysis graphs. You will stay on top of things at all times – both during your experiments and months or years later. Your ZEN Connect projects keep all associated datasets together. It's never been easier to share results and collaborate with others as a team.

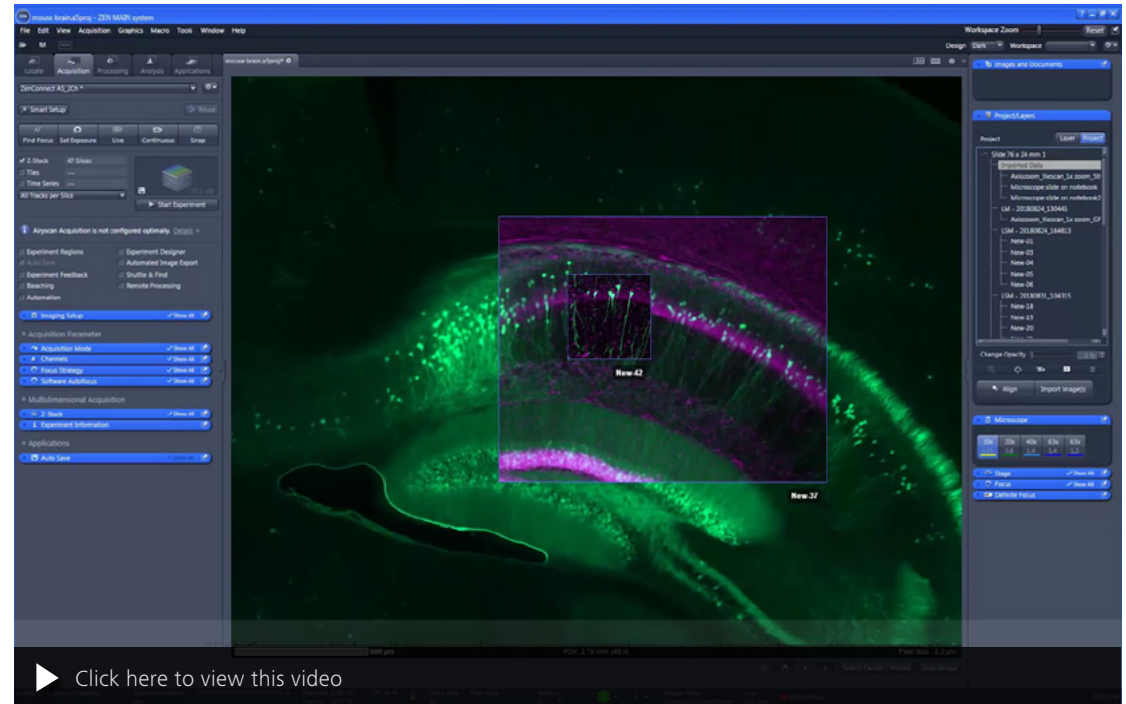

Connect all your imagery: With ZEN Connect you bring images and data from any system or modality together. You always keep the context and the overview about all data from your sample.

The powerful integrated 3Dx1 Viewer, powered by arivis®, is optimized to render the large 3D and 4D image data you have acquired with the fast LSM 980. You can create impressive renderings

and movies for meetings and conferences. After all, a good picture can say more than a thousand words.

Expand Your Possibilities


- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

Get More Data from Your Sample

The real value of microscopy images is in the data they provide. The CZI file format of ZEN microscopy software makes sure that all important metadata of your experiments are safely stored and can be accessed openly for cross-platform data exchange. ZEN provides numerous analysis tools to extract all kinds of information from your images.

Building analysis workflows that adapt to specific applications is not an easy task. It requires knowledge of image processing and the ability to assemble a series of image operations. ZEN addresses this challenge with the BioApps modules for efficient image analysis. Each module is optimized for one type of application, e.g., cell counting or confluency measurement, with tailored segmentation settings and streamlined data presentation. If your applications require customized workflows, the wizard-based ZEN Image Analysis module will guide you step by step to create your unique measurements.


Within an image analysis workflow, segmentation and object classification are two of the most challenging steps. ZEN Intellesis uses the latest machine learning algorithms to make these steps easier and more accurate, also allowing you to execute training on your own data sets. You can integrate the individual models seamlessly into your ZEN image analysis workflow.

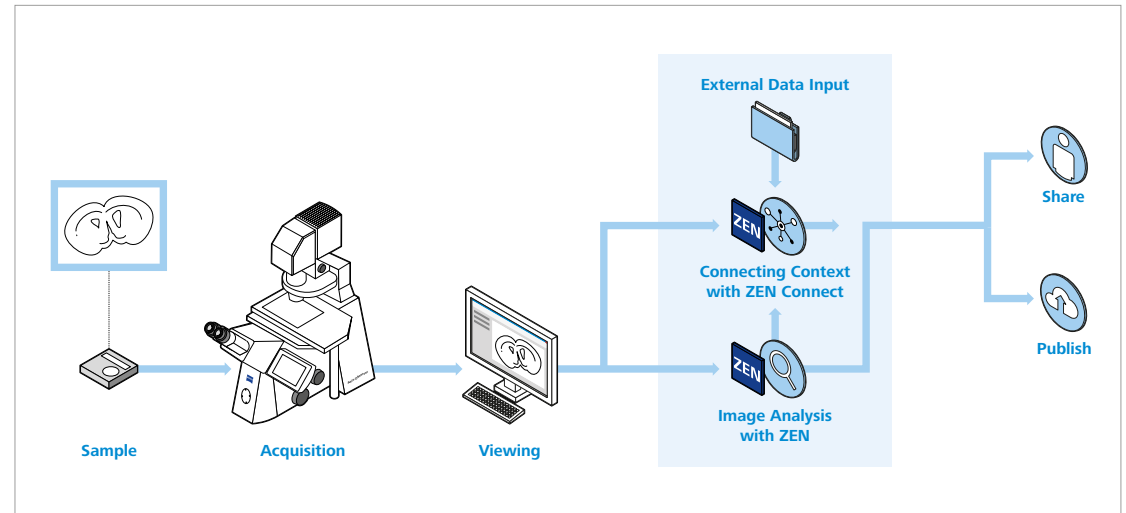

ZEN microscopy software integrates all steps from your sample to reproducible data for publication.



ZEN Intellesis: Use the power of machine learning to easily segment your images.

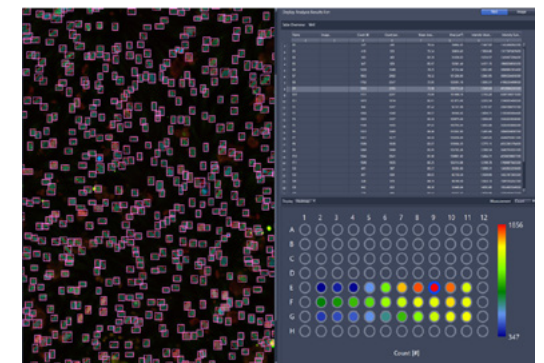

ZEN BioApps: From beautiful images to valuable data – analyze your images efficiently.



Expand Your Possibilities

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service


Combine Multiple Super-resolution Techniques

Combine your LSM 980 with Elyra 7 and Lattice SIM² and you can always choose the best super-resolution technique for your experiment at hand. Lattice SIM² technology brings structured illumination microscopy (SIM) to a new level. Groundbreaking light efficiency gives you gentle live cell super-resolution imaging down to 60 nm resolution with incredibly high speed. With 6 different imaging modes to choose from, even lossless super-resolution imaging becomes possible, meaning you can now get one super-resolution image per acquired raw image for even gentler imaging of your most delicate samples. Add single molecule localization microscopy (SMLM) for techniques such as PALM, dSTORM and PAINT. You can now choose freely among your labels when imaging with resolutions down to 20 nm laterally and 50 nm axially. High power laser lines allow you to image your sample with ease, from green to far red.

Whether in an imaging facility or a single lab, your microscope users will appreciate the wealth of techniques for gentle 3D live cell imaging with super-resolution in one single system.

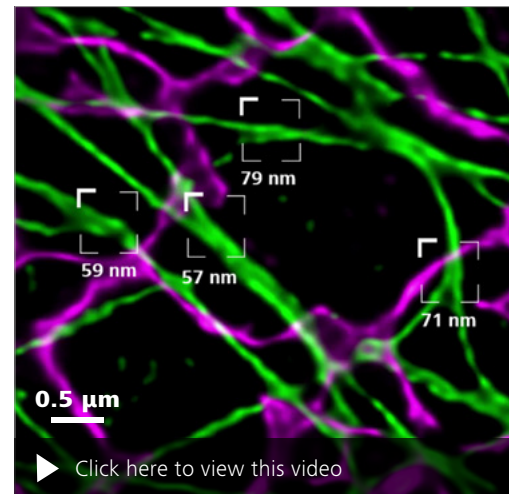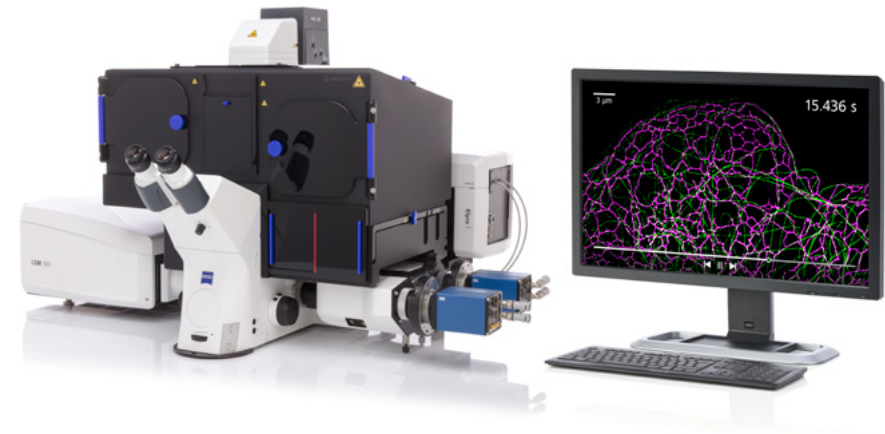

Lattice SIM²: Simultaneous imaging of the endoplasmic reticulum (Calreticulin-tdTomato, magenta) and microtubules (EMTB-3xGFP, green) in a Cos-7 cell reveals highly dynamic interaction of these organelles at resolutions well below 100 nm. Objective: Plan-Apochromat 63x / 1.4 Oil

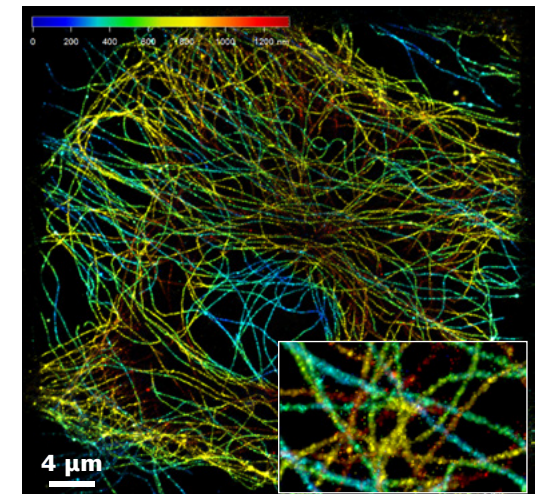

SMLM: With ZEISS Elyra 7 you can image a z-depth of 1.4 µm in a single acquisition. 3D SMLM image of Alexa 647 α-tubulin color coded for depth. Sample courtesy of M. W. Davidson, Florida State University, USA.

Expand Your Possibilities

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

Multiphoton microscopy (commonly referred to as two-photon microscopy or non-linear optical (NLO) microscopy) is well established as the preferred method for non-invasive and deep tissue imaging of living (in vivo and in vitro) or fixed samples and as such has revolutionized in particular neuroscience imaging. Multiphoton microscopy capitalizes on one of the fundamental characteristics of light, that longer wavelengths (600 – 1300 nm) are less absorbed and less scattered by tissues; travelling deeper into the sample while still forming a focal point. The required energy to excite a fluorescent dye is now provided not by one photon but by two photons with half the energy each. And consequently, the probability of two photons to reach the fluorophore at the same time is only high at the focal point.

So all emission light originates from the focal plane and can be efficiently directed to a detector (non-descanned detector, NDD), generating an optical section while omitting a pinhole. Sensitive GaAsP detectors can be placed directly on the microscope in the transmission and reflection beam path, or even directly behind the objective lens, for best signal transfer.

An LSM that shares confocal and multiphoton capabilities allows you to access both technologies, depending on your experiments. You can



Energy diagram of two-photon microscopy


use multiphoton excitation while using the confocal detectors, of course with a fully open pinhole. Make use of the advantages of Airyscan 2 including its Multiplex modes and combine deep tissue penetration with enhanced sensitivity, resolution and speed. Factors that are beneficial for your functional imaging experiments, large volume imaging and screening applications. The 34 Channel Quasar detector for simultaneous spectral scans (lambda scans) further expands the opportunities for your successful multiphoton imaging. Several fluorophores can be excited with one NLO wave-

length, efficiently reducing the light exposure to the sample while speeding up the acquisition process. Linear Unmixing easily allows a clear separation of all collected emission signals.


Even non-stained structures can be visualized with multiphoton excitation by second or third harmonic generation (SHG, THG). SHG effects occur on non-centrosymmetric molecules with predominantly periodic alignment, for example, in striated muscle and collagen.

Expand Your Possibilities

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

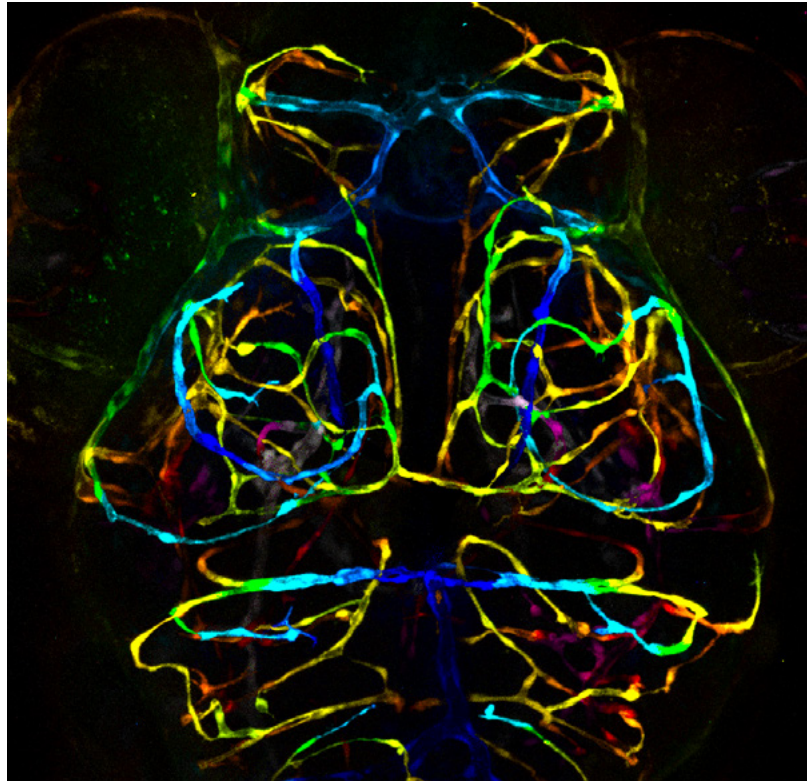

Zebrafish brain vasculature imaged in coronal orientation. Acquired with the two-photon laser excitation at 1,000 nm. The emission light was captured with the GaAsP BiG.2 non-descanned detector. Colour coded 293 µm z-stack. Sample courtesy of the Fish Facility, Leibniz-Institut für Altersforschung – Fritz-Lipmann-Institut e.V. (FLI), Jena, Germany

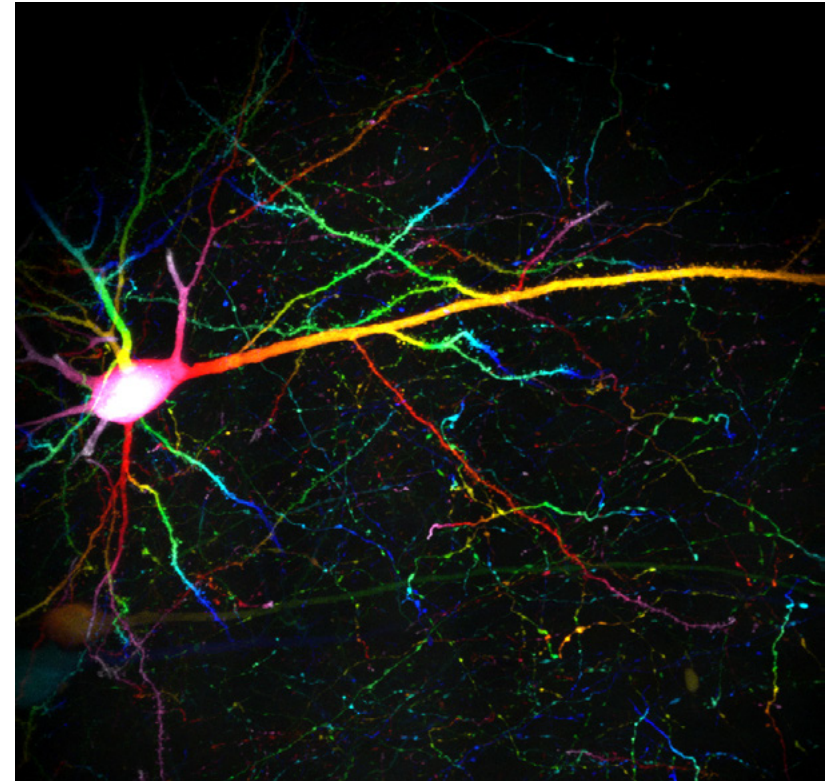

Mouse brain slice with neuronal cytoplasmic GFP label. The 100 µm volume was acquired with two-photon laser excitation at 1,000 nm with the GaAsP BiG.2 non-descanned detector. The dataset was colour coded for depth and an orthogonal projection was created with ZEN blue. Sample courtesy of Prof. J. Herms, LMU, Munich, Germany

Expand Your Possibilities


- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

Image Large Cleared Samples

Tissue clearing opens up a new dimension of optical penetration depth into biological samples such as tissue sections, mouse brains, embryos, organs, spheroids or biopsies.

With Axio Examiner and special objectives—for example, Clr Plan-Apochromat 10×/0.5 nd=1.38, Clr Plan-Apochromat 20×/1.0 Corr nd=1.38 or Clr Plan-Neofluar 20×/1.0 Corr nd=1.45—you can look deep into tissue that has been treated with clearing agents such as Focus Clear or Scale. The cleared tissue becomes almost transparent and the objectives provide the matching refractive index to the immersion medium, delivering crisp contrast. You can now image up to six times deeper than with a multiphoton microscope and up to 60 times deeper than with a conventional laser scanning microscope on uncleared samples.

Get ready to be impressed by the quality of structural information you will retrieve from the deepest layers: expect a big push forward.

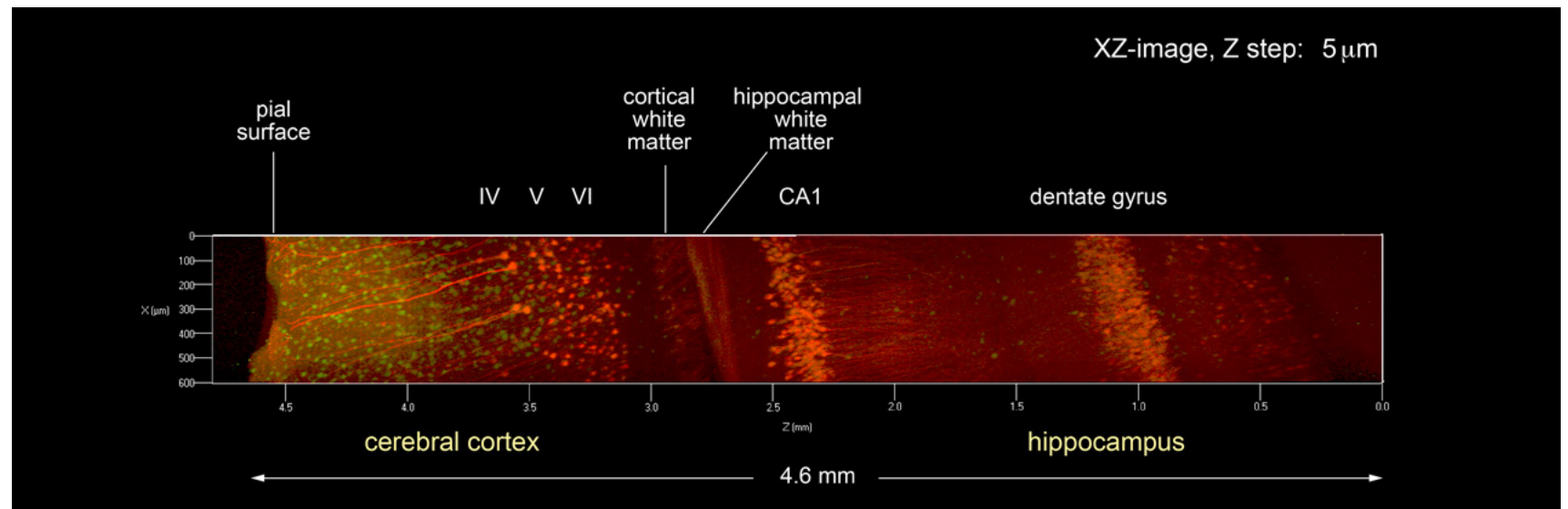

Maximum intensity projection, brain of 7-week old YFP-H mouse, fixed and cleared with Scale clearing technique (Hama et al, Nat Neurosci. 2011).
Courtesy of H. Hama, F. Ishidate, A. Miyawaki, RIKEN BSI, Wako, Japan

Expand Your Possibilities

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service


Correlative Cryo Microscopy: Image the Near-to-native State

Spindle pole bodies are difficult to localize within yeast cells. They are small and rarely occurring structures. ZEISS Correlative Cryo Workflow lets you precisely identify and image such cellular structures in the near-to-native state.

The LSM with the Airyscan detector makes the identification of these structures even easier so further details can be imaged. All images – from a large overview of the entire cell to high-resolution images of these tiny structures – are organized in a ZEN Connect project, providing all data needed to re-locate these cellular structures in the FIB-SEM.




Using ZEISS Crossbeam, TEM lamella of the identified regions can be prepared for cryo electron tomography. Volume imaging is possible as well. Furthermore, the workflow solution allows you to reconnect all data after image acquisition. Images from the Crossbeam or tomograms from the TEM can be combined with the LSM data and can be rendered in three-dimensional context.

Learn more about ZEISS Correlative Cryo Workflow:
www.zeiss.com/cryo

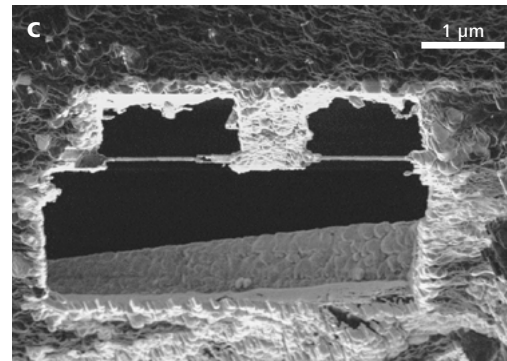

Yeast cells labeled with NUP (nuclear pore complex)-GFP and CNM67-tdTomato. Sample and tomogram courtesy of M. Pilhofer, ETH Zürich, Switzerland

- a) ZEN Connect movie shows the overlay of an LM and EM dataset – from the grid overview to the region of interest identified for further TEM tomography.
- b) Early state of the milling process: Lamella is prepared around the marked region which was identified at the LSM.
- c) FIB image of the prepared lamella; lamella thickness: 230 nm
- d) 3D overlay of the reconstructed and segmented tomogram with LSM dataset (Spindle pole body is false-colored in cyan); nuclear membrane and microtubules were segmented using IMOD.
- e) Segmented and reconstructed tomogram

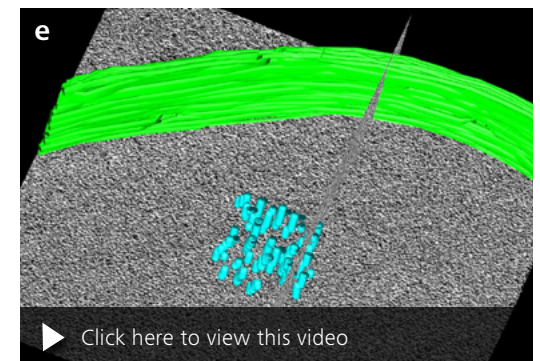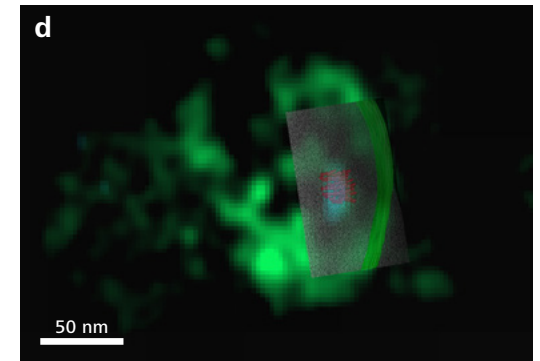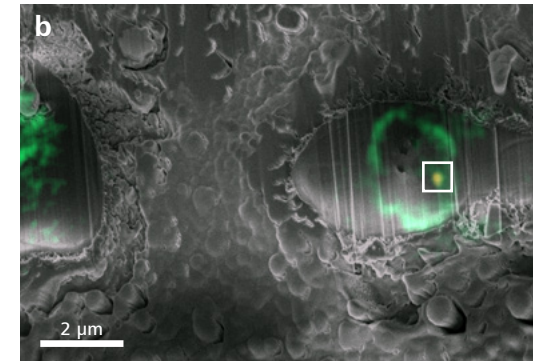


Expand Your Possibilities

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

As your needs grow, your LSM 980 grows with you, forming the basis for a number of enhancements. Like every system from ZEISS, open interfaces and a modular architecture guarantee the seamless interaction of all components now and in the future. These include:



Combine your ZEISS Axio Observer 7 with integrated incubation modules to create the perfect environment for long-term live cell imaging with stable temperature conditions.

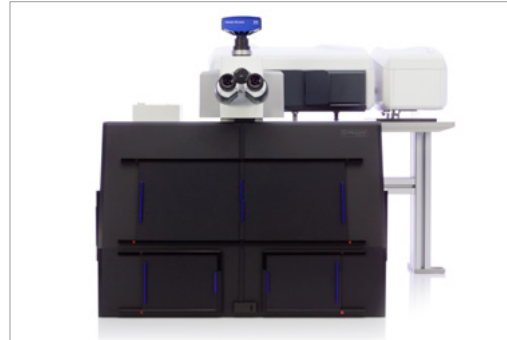

The upright fixed stage microscope ZEISS Axio Examiner.Z1 gives you ample specimen space and room for imaging of whole animals. This stable stand is ideally suited for your demanding multiphoton experiments with incubation for living specimens.




AI Sample Finder automatically detects your sample carrier, adjusts the focus, and finds your sample regions on the coverslip. Even with low-contrast samples, you will access relevant regions with just a click and start your experiment right away.



Enhance your microscope with ZEISS Colibri 7. This flexible and efficient LED light source allows to screen and image your delicate fluorescent samples very gently. You profit from stable illumination and extremely long lamp life.



Add the BiG.2 module with its two GaAsP detectors for FCS, photon counting experiments and FLIM to your ZEISS LSM 980.*

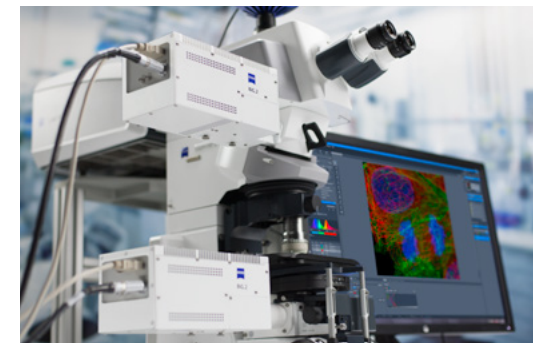


Your BiG.2 works perfectly as a non-descanned detector, also providing a highly sensitive direct coupled detector for FLIM.*

** available upon request*

Expand Your Possibilities


- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service

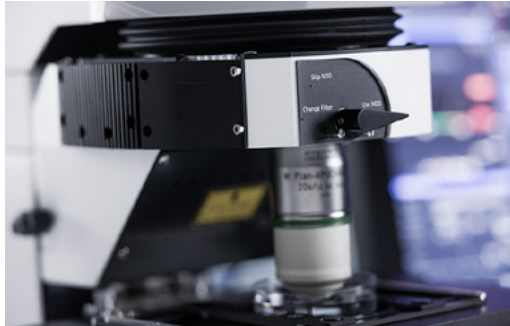


The module GaAsP NDD 2 channels with flexible filter settings completes the ensemble of non-descanned detectors for ZEISS Axio Examiner.Z1.



With autocorr objectives and ZEN microscopy software it's easy to adjust your microscope optics to your sample. You get crisp contrast and better signal to noise – even in your most challenging samples.

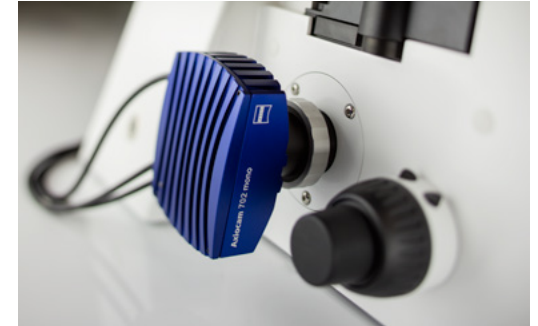


You can add a choice of sensitive ZEISS Axiocams to your ZEISS LSM 980. It's very easy to acquire overview images for your multiposition experiments or to perform light efficient widefield imaging.

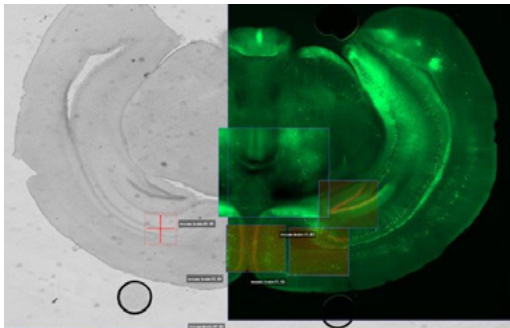

ZEN Connect 2D and 3D Add-on is your gateway to correlative light and electron microscopy (CLEM). Combine the specificity of functional fluorescence imaging with ultrastructural information.



Definite Focus 3 compensates Z-drift and stabilizes the focal position of your sample. You can now perform long-term multiposition and tiling experiments that can last for multiple days.

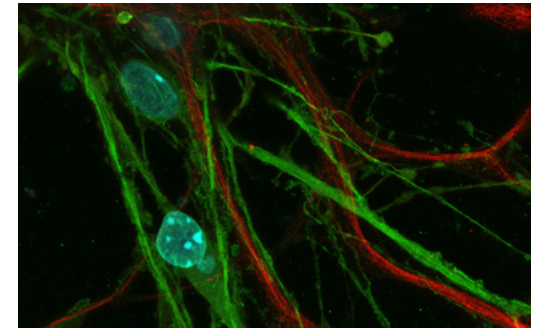

Collect all labels simultaneously with the numerous channels of ZEISS LSM 980 and accelerate the deconvolution process significantly with the CUDA enabled GPUs. Add enhanced resolution and signal to noise to the multi-channel flexibility of imaging with your ZEISS LSM 980.

Expand Your Possibilities

- › In Brief
- › **The Advantages**
- › The Applications
- › The System
- › Technology and Details
- › Service


ZEISS Predictive Service

Maximizes System Uptime

Once connected to your network and activated, this advanced technology will automatically track the health status of your instrument and collect system log files in the background to improve remote diagnosis.

Relevant technical data such as operating hours, cycle counts or voltages are periodically monitored via a secure connection to our data center. The ZEISS Predictive Service application evaluates the performance of your microscope as system data can be received and analyzed.

Our support engineers will diagnose any issues by analyzing data on the Enterprise Server – remotely and without interruption to your operation.

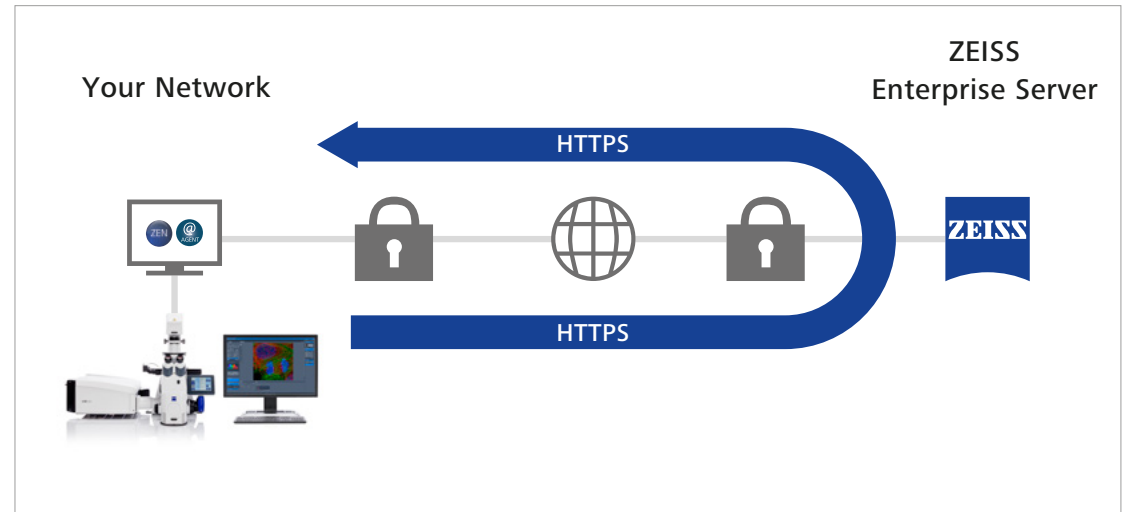

■ Maintain highest system availability

Increase your uptime through close monitoring of the system's condition as remote support can often provide immediate solutions.

■ Data security

Ensure highest data security standards using well established technologies like PTC Thingworx and Microsoft Azure Cloud. No personal or image data is uploaded, only machine data.

■ Fast and competent support

Use secure remote desktop sharing to easily get an expert connected.

■ Optimum instrument performance

As the status of your system is monitored, necessary actions can be planned before they become urgent.

Tailored Precisely to Your Applications

- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details
- › Service

Typical Applications, Typical Samples	Task	ZEISS LSM 980 Offers
Antibody labelled Specimen	Resolve morphological structures down to 120 nm (XY) / 350 nm (Z) at 488 nm excitation while avoiding photobleaching	Airyscan 2 with SR and Multiplex mode for efficient super-resolution imaging
	Acquire large fields of view and tiling experiments	
	Increase your fluorescent range by using labels up into the Near Infrared emission range	Extension of the Quasar detectors with a special dual channel NIR Detector for channel or lambda imaging
Cleared tissue	Investigate structures of your cleared tissue down to 5.6 mm in Z	Special objective corrected for immersion medium of refractive index 1.38 or 1.45 working with confocal or multiphoton imaging on Axio Examiner
Live cell imaging	Study the motility of vesicles and organelles	Airyscan 2 in Multiplex mode for gentle imaging with high frame rates
	Follow fast processes such as Calcium waves, muscle contractions, blood flow, cilia beating while keeping structural information	Airyscan 2 with Multiplex mode for gentle imaging at very high frame rates at confocal resolution
	Acquire information about dynamics and concentration of molecules in living cells	FCS measurement of fluorescent molecules using up to 7 measurement channels.
	Collect additional information of non labelled subcellular structures generating second and third harmonic signals	LSM 980 NLO with broad excitation range up to 1300 nm
Live cell culture with two labels	Study the motility of subcellular structures	Airyscan 2 with GaAsP detector to image 2 colors with time lapse imaging in 2D or 3D at 2.4 frames per second and up to 23 frames per second in Multiplex mode
	Explore the interaction of two proteins with fluorescent lifetime microscopy	BiG.2 as detector for FLIM and third party electronics and software*
	Explore the interaction of two proteins using the Förster Resonance Energy Transfer effect	FRET analysis tool
	Study dynamic characteristics of two or more proteins of interest	FCCS measurement of two or more different fluorescent molecules using up to 7 measurement channels
Live cells with multiple labels	Image over long time in an automated way	Experiment Designer to acquire complex experiments. Combination of different acquisition modes of the LSM system, e.g. spectral imaging, tiling with Airyscan 2 Multiplex mode at super-resolution. Combine all findings with ZEN Connect.
	Reveal mechanistic details in live cells, e.g. moving organelles, vesicle trafficking, membrane reorganization.	Elyra 7 with SIM ²
	Resolve structural details down to 60 nm in 3D and multiple colors.	
Live or fixed samples with multiple labels and overlapping emission signals	Use the fluorescent labels that you need to perform your multi-color experiment to investigate the interplay of multiple proteins.	Parallel acquisition of all signals with spectral imaging at 5 frames per second and online or post processed linear unmixing
Cellular structures with weak labels	Image subcellular structures at physiological expression levels	LSM 980 with Airyscan 2 with GaAsP detectors

(* available upon request)

Tailored Precisely to Your Applications

- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details
- › Service



Typical Applications, Typical Samples	Task	ZEISS LSM 980 Offers
Living organisms/animals	See the interaction of cells within living tissue	Multiphoton extension of LSM 980 to LSM 980 NLO
	Imaging of living tissue with cells expressing multiple different fluorescent proteins	LSM 980 NLO with two NLO excitation lines and up to 7 non descanned channels in reflection for spectral imaging
	Follow fast processes such as Calcium waves, muscle contractions or blood flow, in whole organisms or tissue culture while keeping structural information.	Combine multiphoton excitation with Airyscan 2 Multiplex to capture fast processes in deep tissue while maintaining structural information and avoid phototoxicity
	Volume defined photomanipulation of cellular and subcellular structures	LSM 980 NLO with high power and broad wavelength range for efficient dye manipulation
Plant roots	Follow the changes of subcellular structures over time with high resolution	Airyscan 2 with GaAsP detector for super-resolution imaging beyond 40 µm deep into tissue with up to 47 frames per second (512×512)
Model organisms, e.g. Zebrafish, <i>Drosophila</i> or <i>C. elegans</i>	See fine details of the organization and dynamics of endogenously expressed fluorescent proteins	Airyscan 2 with GaAsP detector for super-resolution imaging beyond 40 µm deep into tissue
	Image large fields of view at high volume rate to capture developmental processes	Airyscan 2 with Multiplex mode for high frame rates at confocal resolution
	Capture 3D volume image of total specimen	LSM 980 NLO with two NLO excitation lines for volume imaging of multiple labelled samples

ZEISS LSM 980 at Work

- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details
- › Service

Expand the number of labels

To capture the complex world of biology, the ability to expand the number of labels is a great advantage. LSM 980 can image simultaneously multiple labels, covering a wide emission range up to 900 nm. These Cos-7 cells were labelled with 4 different fluorophores, two of which have their emission peak in the near infrared range (NIR), Alexa 700 and Alexa 750. Utilizing the flexible LSM 980 Quasar and NIR detectors, all labels were imaged with optimal sensitivity.

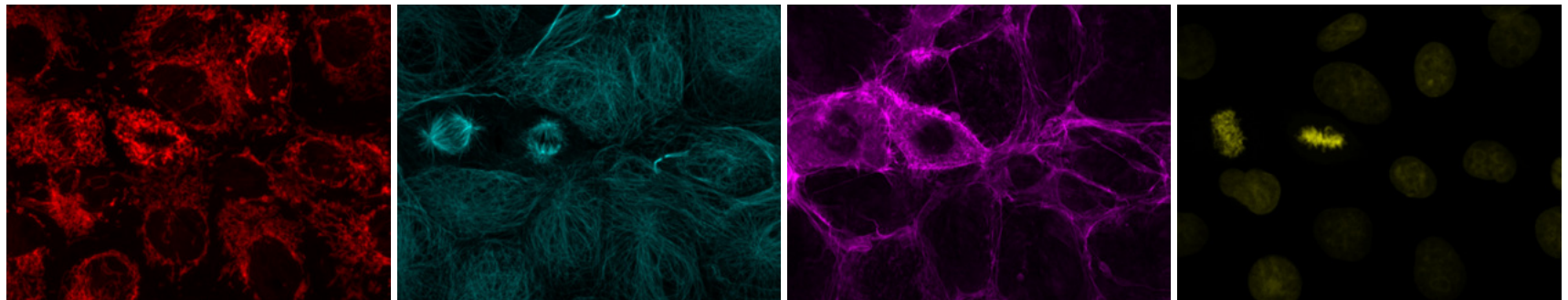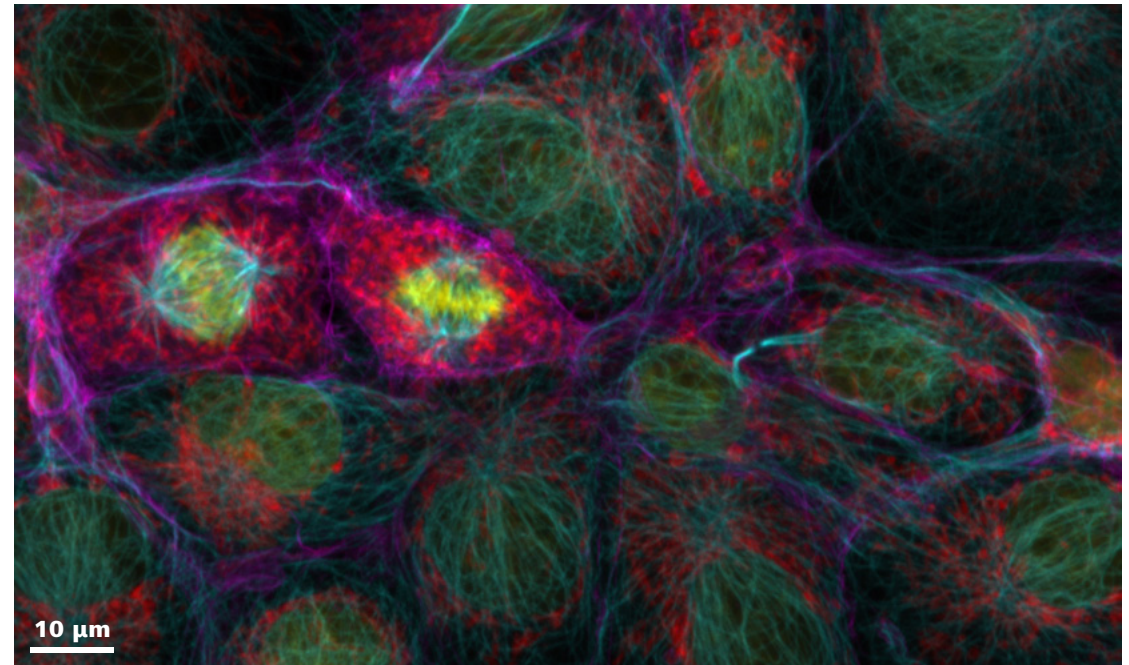

Cos-7 cells Anti-TOM20 AF750 (red), Anti-Tubulin AF568 (cyan), Actin Phalloidin-OG488 (magenta), DAPI (yellow). Imaged with LSM 980 including the ZEISS NIR detector in channel mode. The fluorescent signals were separated by Linear Unmixing, facilitating clear separation between the spectrally overlapping dyes Alexa 700 and Alexa 750. Sample courtesy of Urs Ziegler and Jana Doehner, University of Zurich, ZMB, Switzerland.


ZEISS LSM 980 at Work

- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details
- › Service

Navigate and correlate with ease

As the world of microscopy transitions gradually to larger samples, it becomes more important to maintain the positional context and keep a record of the areas captured. AI Sample Finder automatically classifies the sample carrier, identifies the sample, finds the focus, and creates a fast overview image using the T-PMT detector or camera. You can freely navigate using the overview image for orientation, and effortlessly move to the structures of interest. Making sure you only spend time imaging regions that hold information for your research. ZEN Connect correlates all data associated with the sample.


In this example, mouse intestinal tissue was labelled with three fluorophores covering an emission spectrum of 500 – 850 nm. AI Sample Finder automatically identified the carrier and created an overview image using the T-PMT to capture the Alexa 488 label. The overview image is used for sample navigation and identification of regions of interest. The ZEISS LSM 980 Quasar and NIR detectors were used to acquire images of the visible and invisible labels with optimal sensitivity.

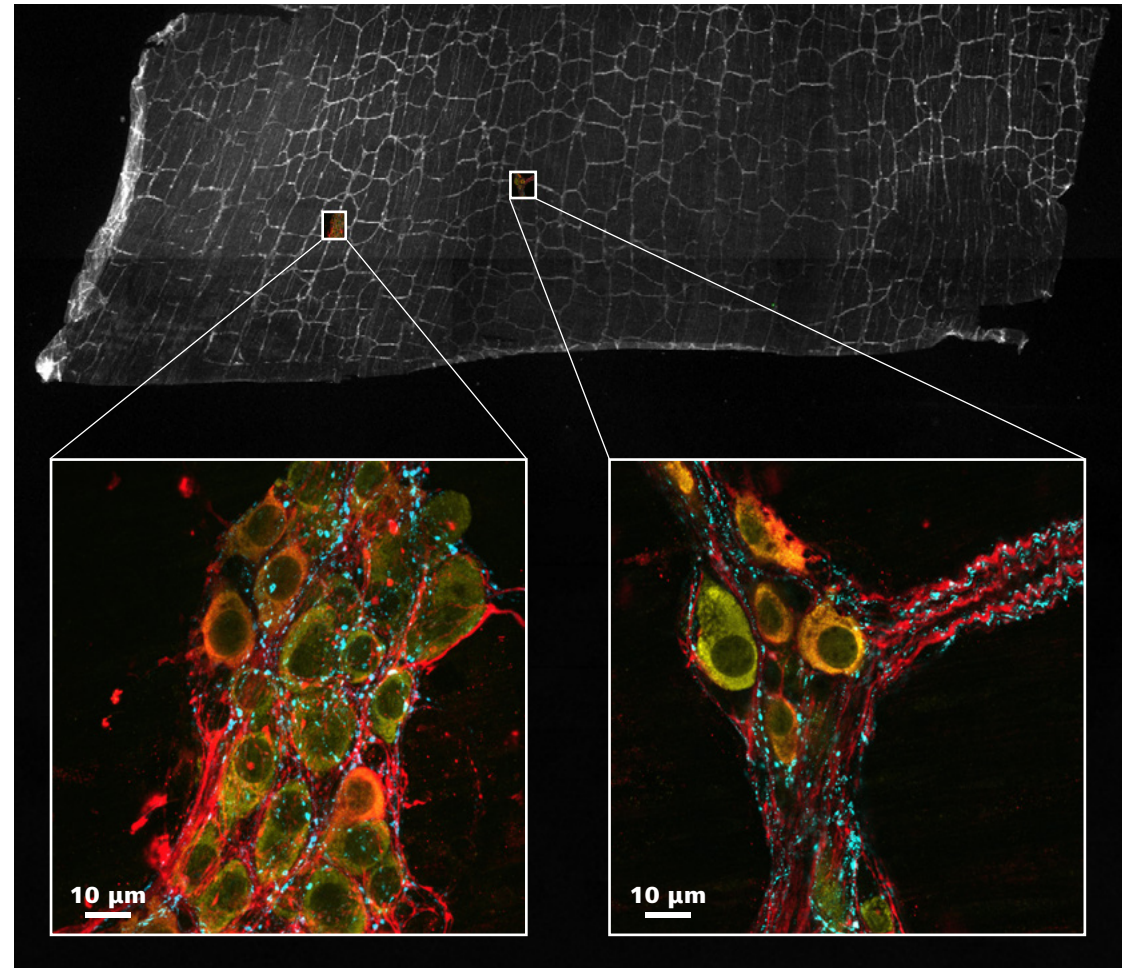

Mouse intestine tissue section stained for Substance P (cyan, Alexa 488) labeling the presynaptic contacts in the enteric nervous system, HuC/D (yellow, Alexa 568) labeling the enteric neurons, and neuronal Nitric Oxide Synthase (nNOS, red, Alexa 750) labeling a sub-population of enteric neurons. Sample Courtesy of Pieter Vanden Berghe, LENS & CIC, University of Leuven, Belgium.

ZEISS LSM 980 at Work

- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details
- › Service

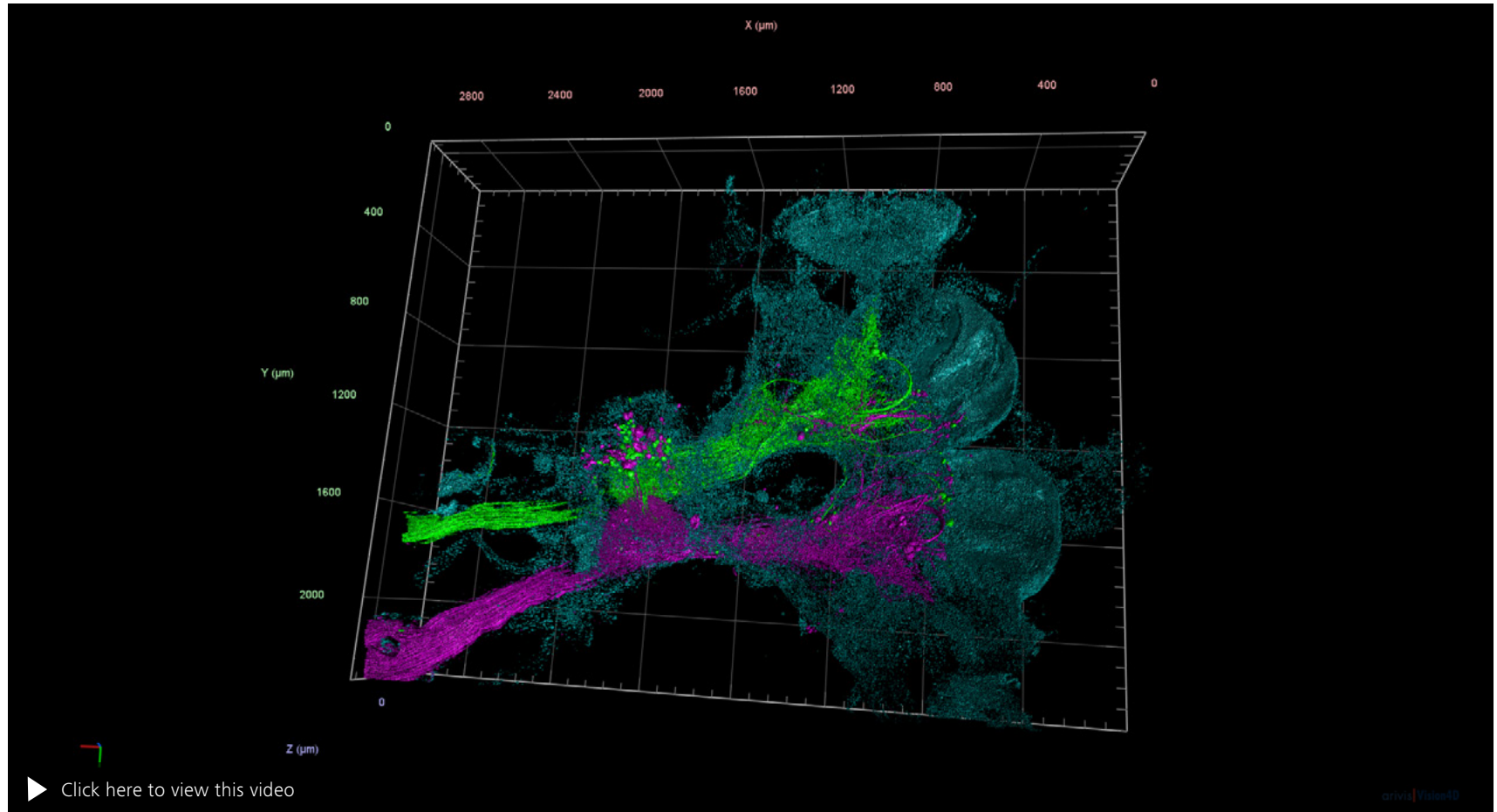

The brain, thoracic and abdominal ganglia of the cockroach are joined together by bilateral connective bundles of ascending and descending interneurons forming the ventral nerve cord. In this preparation, left and right connectives were individually labelled (Alexa 488: green, Alexa 647: magenta) posteriorly to the suboesophageal ganglion to observe the extension of their innervation within the different neurophils, and throughout the ipsi- and contralateral parts of the brain (DNA labelled with DAPI: cyan). Imaging was performed using Tiling and Stitching to capture the complete volume (3×2.3× 0.26 mm). 3D animation of the complete dataset was done with arivis Vision 4D, ideal for rendering and analyzing large datasets. The 4D viewer in arivis Vision 4D can be configured to adjust the appearance of individual channels independently to highlight specific features. These settings, along with clipping planes or the varying opacity of individual channels, can be stored into key frames which the software automatically interpolates between to produce a seamless animation.

These animations can be previewed and edited prior to producing high resolution video renders. Sample courtesy of M. Paoli, Galizia Lab, University of Konstanz, Germany

ZEISS LSM 980 at Work

- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details
- › Service

Oocytes store all the nutrients to support early embryonic development, and are therefore very large cells with a large nucleus. Oocytes need to divide before fertilization. How to make cell division work in this very large cell is the topic investigated by P. Lenart's lab.

They have shown that, surprisingly, an actin network is required to collect chromosomes scattered in the oocyte nucleus. They are then handed over to microtubules, which capture chromosomes and align them on the spindle. The actin-driven and microtubule-driven transport phases have very different speeds and show other differentiating characteristics that can be distinguished by tracking chromosome motion.


Peter Lenart says: "This is a nice imaging challenge, because chromosomes are scattered in the spherical nucleus with a diameter of 80 μm and are transported over a period of approximately 15 minutes. Back in 2005 we could acquire stacks every 45 s, which was sufficient to distinguish actin- and microtubule-driven phases. Using the new, high resolution trajectories shown here we hope to learn about the details of the transport mechanism."



Meiosis in starfish oocytes

The depth coding shows a subset of 52 μm . The movie shows the transport of chromosomes, labeled by Histone 1-Alexa 568, in a starfish oocyte undergoing meiosis.

A z-stack of 67 μm was acquired every 2.4 seconds with Airyscan CO-8Y mode. Concomitant with chromosome transport, the nucleolus (the large spherical structure) is disassembling. Courtesy of P. Lenart, MPI for Biophysical Chemistry, Göttingen, Germany.

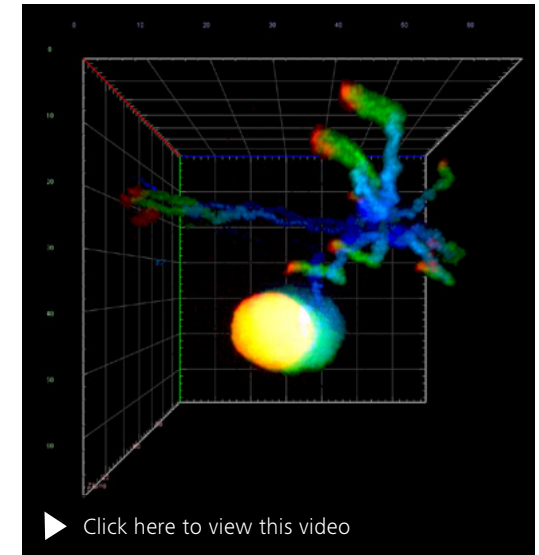

Meiosis in starfish oocytes


The rendering is a projection of the process along z-axis (maximum intensity) and time (color-coded projection); to illustrate the movement of the chromosomes within the volume of the nucleus.

Reference:


- Lenart P, et al. Nature. 2005 Aug 11;436(7052):812-8.
Mori M, et al. Curr Biol. 2011 Apr 12;21(7):606-11.
Bun P, et al. Elife. 2018 Jan 19;7. pii: e31469. doi:10.7554/eLife.31469.
Burdyniuk M, et al. J Cell Biol. 2018 Aug 6;217(8):2661-2674.

ZEISS LSM 980 at Work

- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details
- › Service

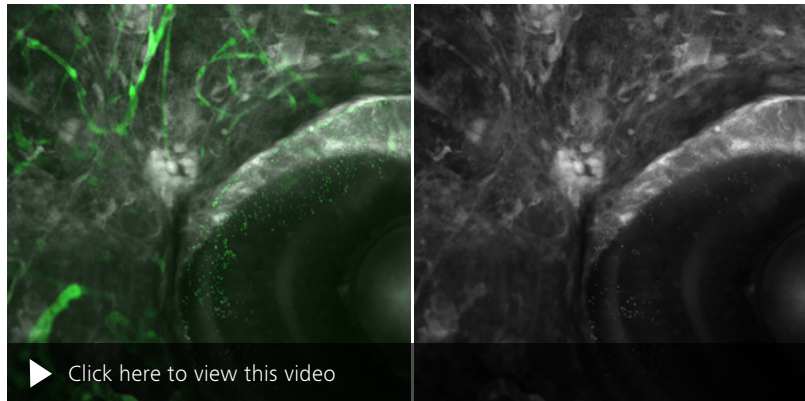

Zebrafish brain and eye vasculature (green) and Second Harmonic Generation (grey) in sagittal orientation. A volume of 267 μm was acquired with the two-photon laser at 1,000 nm and emission was detected with the GaAsP BIG.2 detector. SHG allowed the visualization of the tissue structures, such as the retinal cells and ocular muscles. Sample courtesy of the Fish Facility, Leibniz-Institut für Alternsforschung – Fritz-Lipmann-Institut e.V. (FLI), Jena, Germany

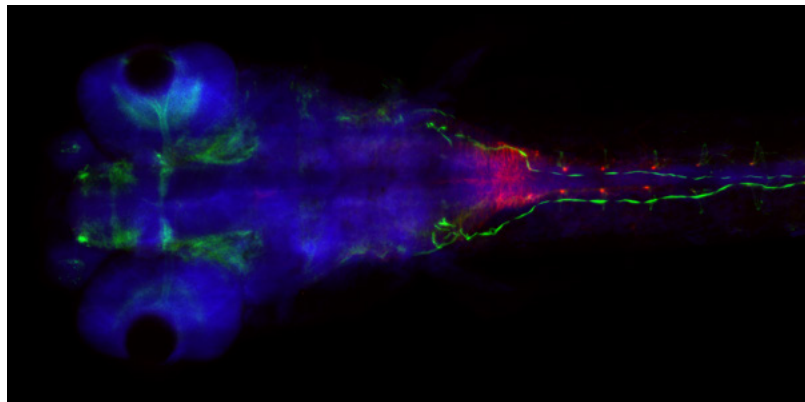




Zebrafish larvae immunolabelled for nervous system (Alexa-488), cell nuclei (DAPI), acetylated tubulin (Alexa-633). All three dyes were excited with the two-photon laser at 780 nm and emission was detected with a lambda stack in 3D Tiling (157 μm , 15 tiles). The fluorophore spectra were separated with Linear Unmixing and the 3D tiles were Stitched in ZEN Blue. This image is the orthogonal projection of the 3D dataset. Sample courtesy of Dr. H. Reuter, Molecular Genetics Group, Leibniz-Institut für Alternsforschung – Fritz-Lipmann-Institut e.V. (FLI), Jena, Germany

Zebrafish is a well established model for studying development of the vascular system. Multiphoton imaging is a great way to capture the intricate vasculature patterns in zebrafish brain at great depth. Additionally, through the Second Harmonic Generation (SHG), structural information of the surrounding tissues can be captured without the need of additional labelling.

To minimize acquisition time and photobleaching it is useful that you can excite multiple fluorophores with a single two-photon laserline, allowing you to image multiple labels simultaneously. With LSM 980 and the 32-Channel Quasar detector you can perform lambda scans, collecting and unmixing the emitted signals over the whole visible spectrum. In this example, all 3 labels (Alexa-488, Alexa-633 and DAPI) were excited and detected simultaneously. 3D Tiling and Stitching was used to capture the large sample volume, maintaining the context in the tissue. LSM 980 offers you the best of both worlds, allowing you to combine the strength of the NLO excitation, with the advantages of spectral detection, and all the imaging modalities of ZEN blue.

ZEISS LSM 980 at Work

- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details
- › Service

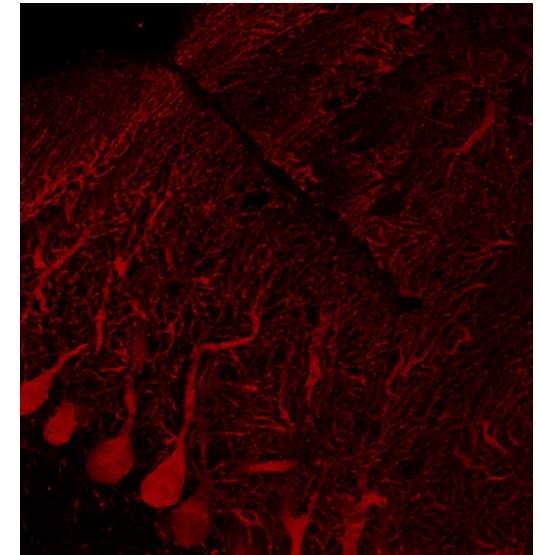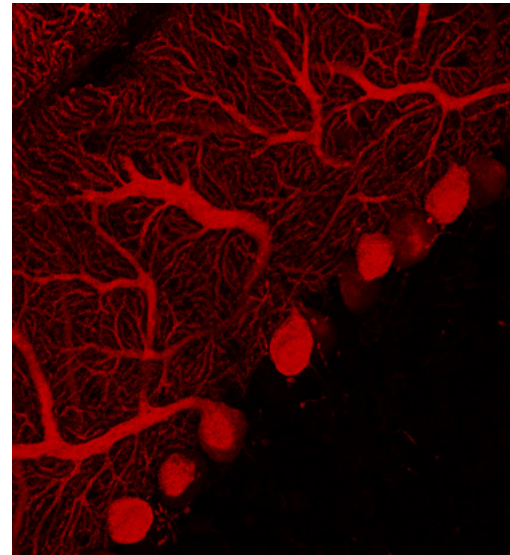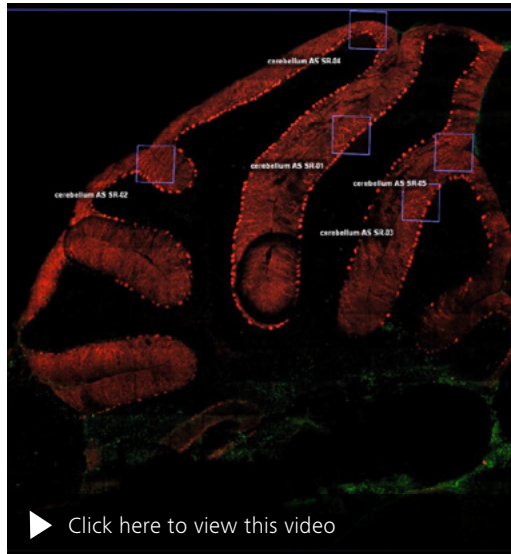

Mouse brain cerebellum labelled with anti-calbindin (Alexa-568) and anti-GFAP (Alexa-488). The fluorophores were both excited with the two-photon laser at 780 nm and the emission spectra were simultaneously collected by the BIG.2 detector. 3D Tiling and Stitching were used to cover whole structure, and an orthogonal projection was created in ZEN Blue. Specific areas of interest were imaged with the Airyscan 2 detector in order to acquire high resolution images of the Purkinje cells. The Airyscan 2 datasets were processed and orthogonal projections were created with ZEN Blue. The individual super-resolution images were aligned with the cerebellum using ZEN Connect. Sample courtesy of L. Cortes, University of Coimbra, Portugal.


Multiphoton microscopy can be combined with 3D Tiling and Stitching in order to image large samples, such as this example of mouse cerebellum. Airyscan 2 imaging in super-resolution mode

can be used to acquire super-resolution images of specific areas of interest, and can be seamlessly combined with two-photon imaging. ZEN Connect can bring all the information from your


different experiments together, allowing you to map the high-resolution images on the larger structure, maintaining the context and simplifying your file organization.

ZEISS LSM 980 at Work


- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details
- › Service

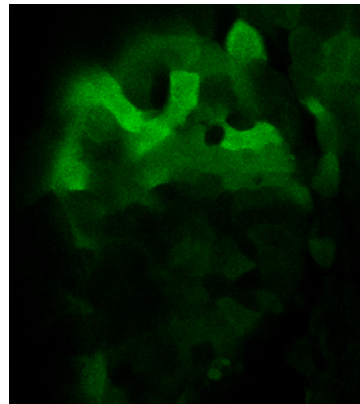


A) YFP intensity (bottom photobleached)

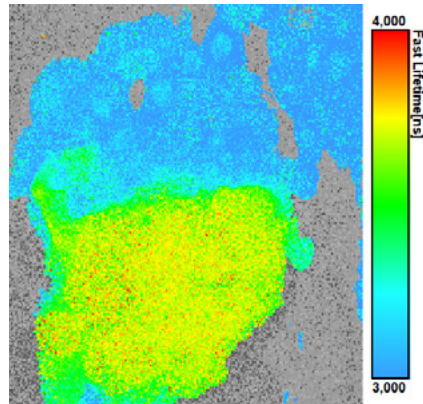

B) Donor / CFP lifetime (bottom photobleached)

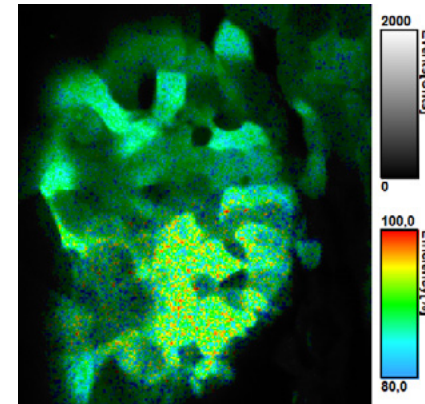

C) FRET efficiency based on Donor lifetime (bottom photobleached)

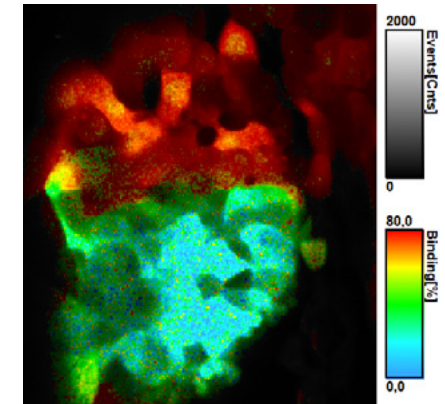

D) Binding between CFP and YFP based on Donor lifetime (bottom photobleached)

Fluorescence Lifetime Imaging (FLIM)

The lab of Marcos Gonzalez-Gaitan is investigating the role of small GTPase during Zebrafish embryonic development. The focus of their work lies in identifying when and where these GTPases are active during the oriented division of ectodermal progenitor (epiblast) cells. This activity can be monitored by using Förster Resonance Energy Transfer (FRET), in which energy transfer from one chromophore (Donor) to another (Acceptor) only occurs when the two chromophores are closer than <10 nm. By measuring the fluorescent lifetime of the Donor (FLIM-FRET), relevant information can be collected.

In this example, the small GTPase Rac protein was fused to variants of CFP and YFP, as an

intramolecular FRET-pair biosensor to monitor GTPase activity. When the acceptor fluorophore is bleached (Fig.A: lower region of the image), the lifetime of the Donor fluorophore is increased in the same region (Fig. B). The FRET Efficiency is not influenced by the bleached Acceptor fluorophore and stays unchanged for the remaining FRET-pairs (Fig. C). Additional information is given in the Binding fraction (Fig. D), which holds quantitative spatial and temporal information of the currently active FRET pairs.

Quantitative FLIM-FRET analysis allows determining the spatial and temporal activity of two or more interacting molecules. In contrast to measuring FRET by photobleaching or intensity ratio


imaging, lifetime imaging allows precise quantification of FRET Efficiency. In addition, FLIM-FRET is used to quantify the binding fraction for a particular intermolecular FRET pair and the fraction of active sensors for an intramolecular FRET biosensor when using a suited FRET pair.

Data was obtained with LSM Systems equipped with PicoQuant FLIM & FCS kit and the PicoQuant FLIM module in ZEN microscopy software. Multi-photon excitation at 840 nm was used for lifetime measurements of CFP. Acceptor photobleaching was performed with 514 nm laser. The analysis was done within PicoQuant's SymPhoTime64.


*available upon request

ZEISS LSM 980 at Work


- › In Brief
- › The Advantages
- › **The Applications**
- › The System
- › Technology and Details
- › Service

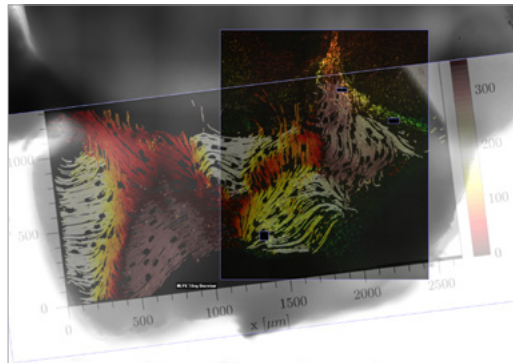

This ZEN Connect project documents the experiment performed with the tissue explant of ependyma from the ventricular system of a mouse brain. All acquired data of the experiment session is kept in context. The overview images by camera and LSM allow to precisely record the localization of the acquired ciliary beating within the sample. The flow map of cilia generated flow along the ependymal wall is added as a reference.

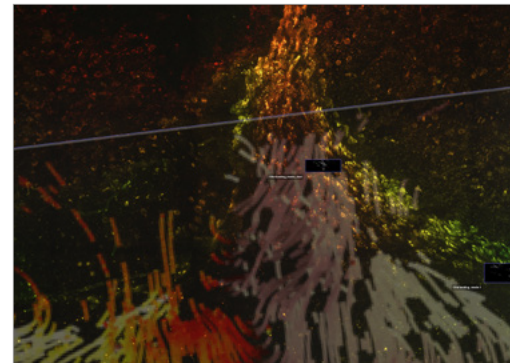

An overview of fluorescently labeled motile cilia on ependyma tissue explant from the mouse brain is quickly acquired by tiling with Airyscan 2 in Multiplex CO-8Y mode to find regions of interest. Z-Stack displayed in colored depth coding. The exact position of the recorded motile cilia is documented.

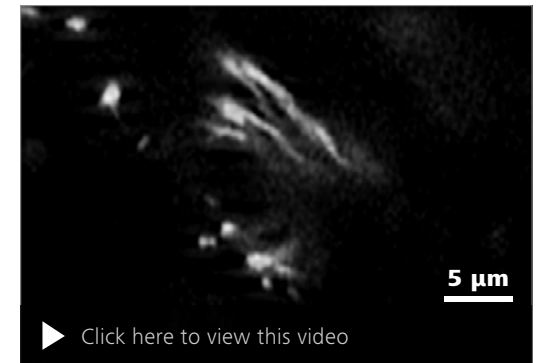

Live imaging with 143 frames per second of fluorescently labeled motile cilia of brain ependyma. Acquired with Airyscan CO-8Y mode combining image quality and speed; for detailed analysis of ciliary beating direction and frequency.

Reference for all images:

G. Eichele, Department of Genes and Behavior, Max Planck Institute for biophysical Chemistry, Göttingen, Germany

Your Flexible Choice of Components

- › In Brief
- › The Advantages
- › The Applications
- › **The System**
- › Technology and Details
- › Service



1 Microscope

- Inverted stand: Axio Observer
- Upright stand: Axio Examiner, Axio Imager
- Port for coupling of Elyra 7
- Camera port
- AI Sample Finder for Axio Observer
- Manual or motorized stages
- Incubation solutions
- Fast Z piezo inserts
- Definite Focus

2 Objectives

- C-Apochromat, C Plan-Apochromat
- Plan-Apochromat
- W Plan-Apochromat, Clr Plan-Apochromat,
- Clr Plan-Neofluar
- LD LCI Plan-Apochromat

3 Illumination

- V laser: 405 nm
- VIS + NIR laser: 445 nm, 488 nm, 514 nm, 543 nm, 561 nm, 594 nm, 639 nm, 730 nm
- Laser for multiphoton imaging: Ti:Sa (single-line laser), InSight X3 / X3+ and Discovery NX (dual-line laser)

4 Detection


- 3, 6, or 34 descanned spectral channels (GaAsP and multialkali PMT)
- NIR Detector (2 channels) with near infrared optimized GaAsP and GaAs detector
- Airyscan 2 detector with optional Multiplex module
- 2 additional GaAsP channels (BiG.2)
- Up to 6 non-descanned GaAsP detectors
- Up to 12 non-descanned GaAsP and multialkali PMT detectors
- Transmitted light detector (T-PMT)

5 Software

- ZEN microscopy software, highlighted modules: Tiles & Positions, Experiment Designer, FRAP, FRET, FCS, RICS, Deconvolution, ZEN Connect, Direct Processing (incl. Denoising and Deblurring) 3Dxl Viewer – powered by arivis®





System Overview

- › In Brief
- › The Advantages
- › The Applications
- › **The System**
- › Technology and Details
- › Service

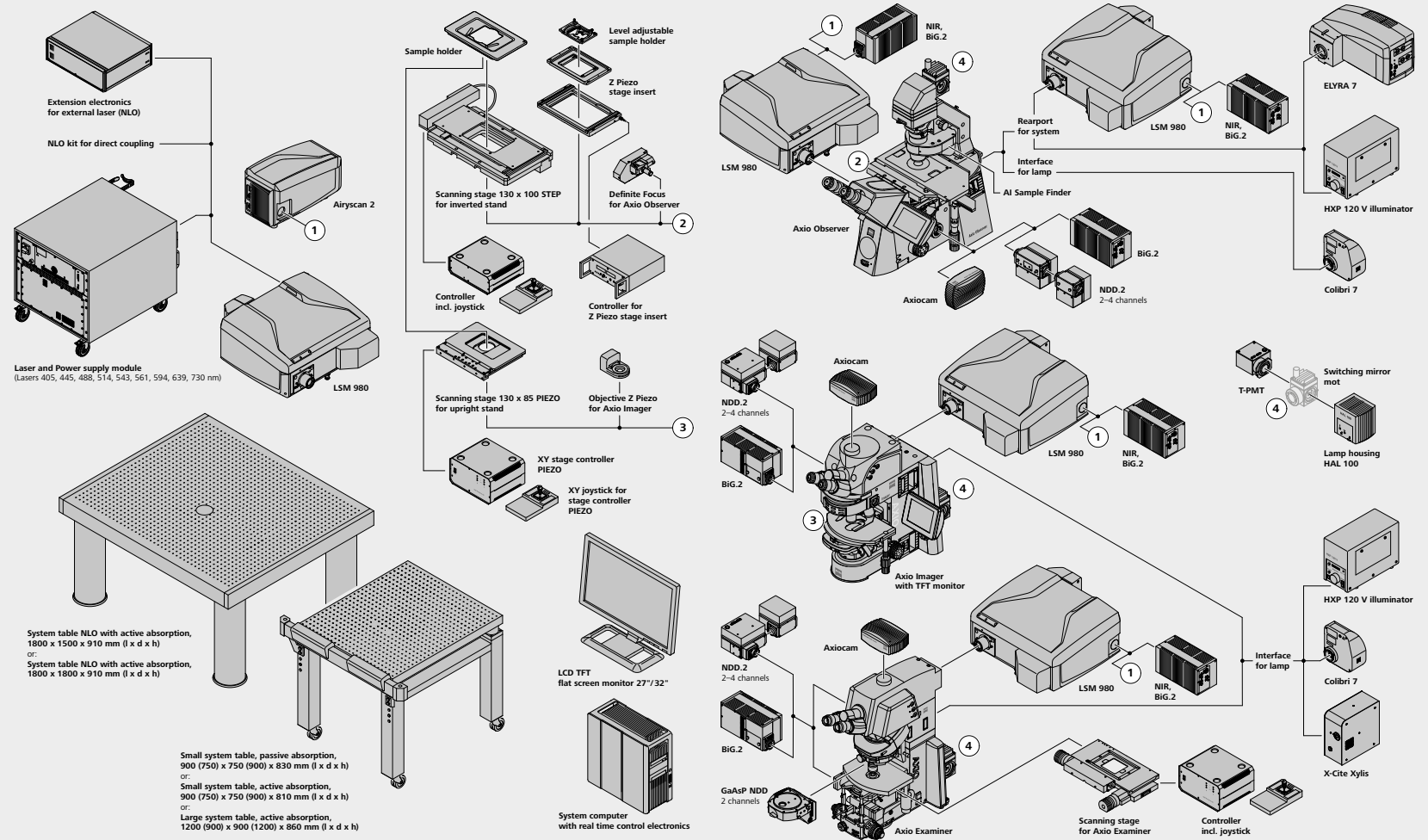

Technical Specifications


- › In Brief
- › The Advantages
- › The Applications
- › **The System**
- › Technology and Details
- › Service

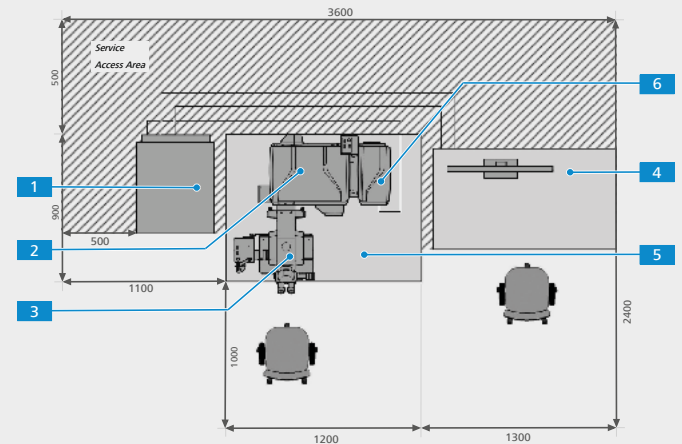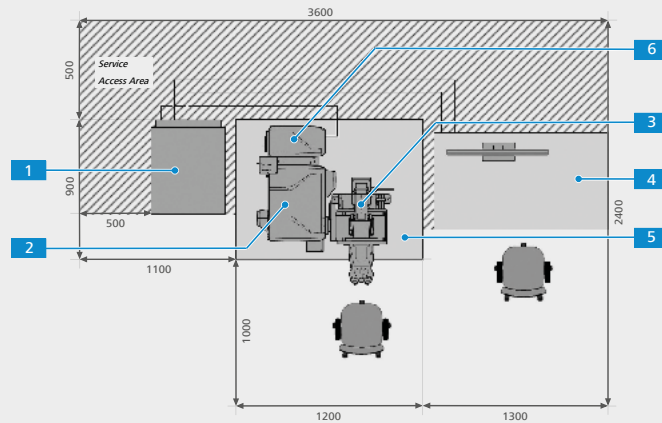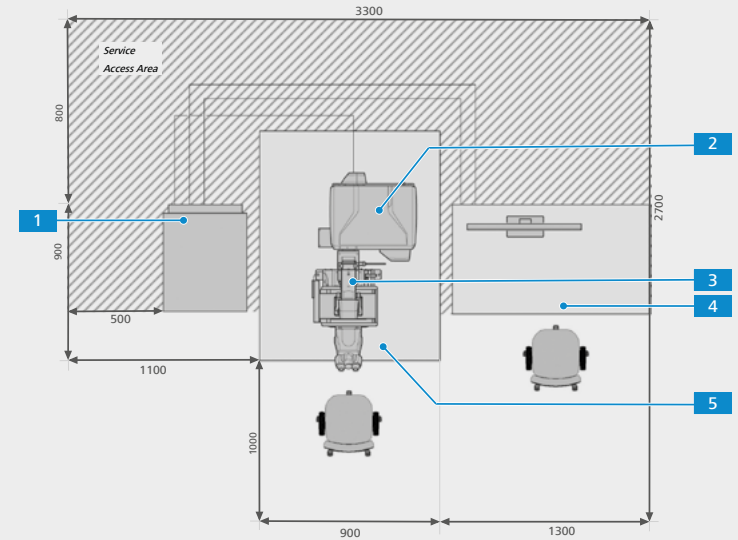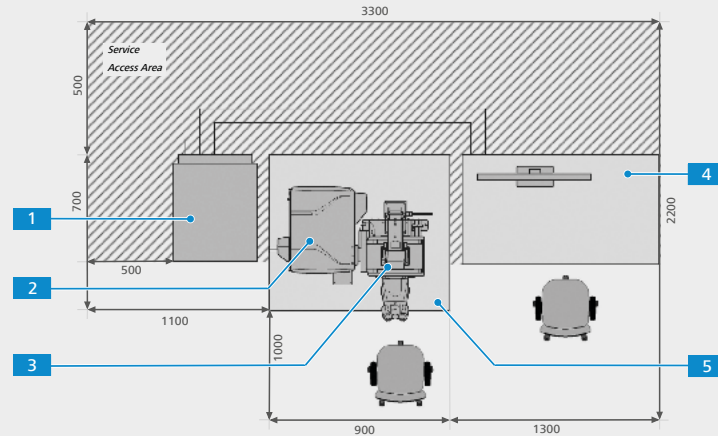


1 Laser and Power Supply Module 2 LSM 980 Scanhead 3 Microscope Stand (Axio Observer, Axio Imager or Axio Examiner) 4 Computer Table 5 System Table 6 Airyscan 2

Technical Specifications


- › In Brief
- › The Advantages
- › The Applications
- › **The System**
- › Technology and Details
- › Service

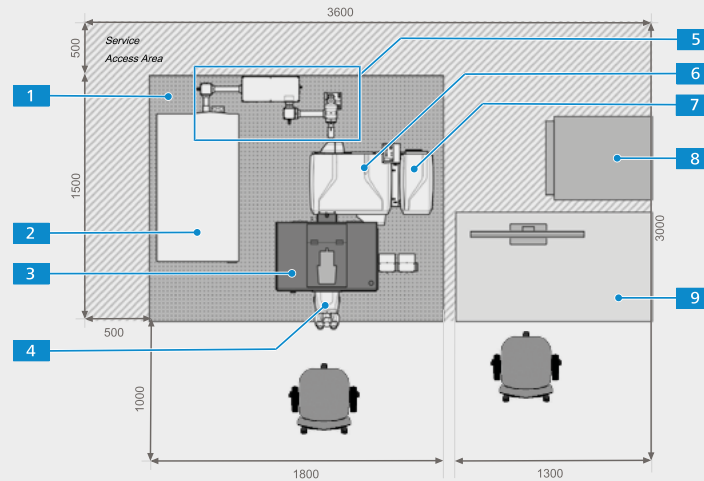

- 1 System Table
- 2 Multiphoton Laser
- 3 Incubation
- 4 Microscope Stand (Axio Observer, Axio Imager or Axio Examiner)
- 5 Laser coupling with AOM for multiphoton laser
- 6 LSM 980 scanning module
- 7 Airyscan 2
- 8 Laser and Power Supply Module
- 9 Computer Table

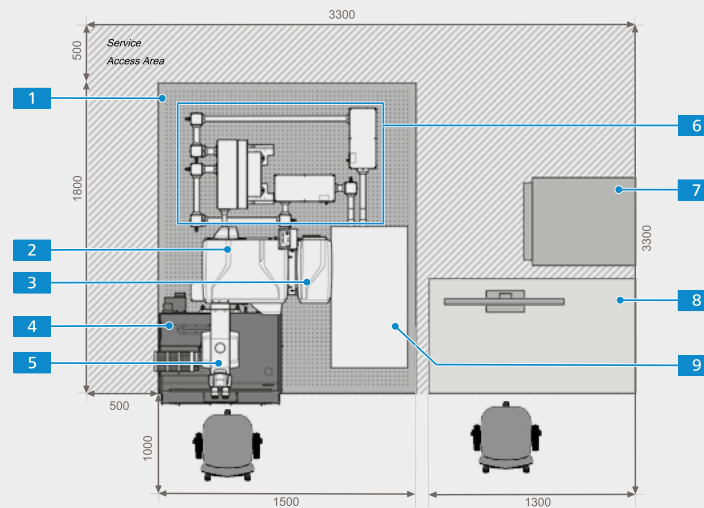

- 1 System Table
- 2 LSM 980 scanning module
- 3 Airyscan 2
- 4 Incubation
- 5 Microscope Stand
- 6 Laser coupling with AOM for multiphoton laser
- 7 Laser and Power supply module
- 8 Computer Table
- 9 Multiphoton Laser

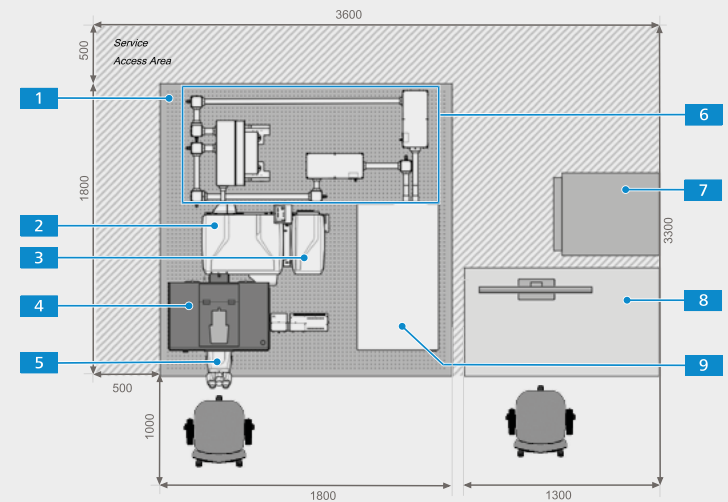

Technical Specifications

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › **Technology and Details**
- › Service

Physical Dimensions	Length (cm)	Width (cm)	Height (cm)	Weight (kg)
Small Passively Damped System Table	90	75	83	130
Small Actively Damped System Table	90	75	81	130
Large Actively Damped System Table	120	90	86	180
Active Anti-Vibration Table (NLO)	180	150	91	475
Active Anti-Vibration Table (NLO)	180	180	91	515
Scanning Module LSM 980	55	45	22	27
Microscope	47–80	29–39	70–72	37–47
Laser and Power Supply module	60	50	56	70
Airyscan 2	40	20	24	12
Fiber Optic Cable, UV	400			
Fiber Optic Cable, VIS	400			
Cables	250			
Microscopes				
Stands	Upright: Axio Imager.Z2, Axio Examiner.Z1 Inverted: Axio Observer 7 with side port or rear port, AI Sample Finder (optional)			
Z Drive	Smallest increment Axio Imager.Z2: 10 nm; Axio Observer 7: 10 nm; Axio Examiner: 25 nm; fast piezo objective or stage focus available; Definite Focus for Axio Observer 7			
XY Stage (optional)	Motorized XY scanning stage, for Mark & Find function (XYZ) as well as Tile Scan (Mosaic Scan); smallest increment of 0.25 µm (Axio Observer 7), 0.2 µm (Axio Imager.Z2) or 0.25 µm (Axio Examiner.Z1)			

Technical Specifications

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › **Technology and Details**
- › Service

Scanning Module

Scanner	Two independent, galvanometric scanning mirrors with ultrashort line and frame flyback
Scanning Resolution	32 × 1 to 8,192 × 8,192 pixels, also for multiple channels, continuously adjustable
Scanning Speed	At 512 × 512 pixels: confocal – up to 13 fps; Airyscan SR – up to 4.7 fps; Multiplex SR-4Y – 25 fps; Multiplex SR-8Y – 47.5 fps; Multiplex CO-8Y – 34.4 fps 19 × 2 speed levels for confocal; 512 × 16 pixels up to 425 fps; up to 6830 lines/sec. 13 × 2 speed levels in Multiplex mode; up to 25 fps for 904 × 904 pixels; up to 17.8 fps at 1,024 × 1,024 pixels
Scanning Zoom	0.6× to 40×; digitally adjustable in increments of 0.1 (Axio Examiner: 0.67 × to 40×)
Scanning Rotation	Can be rotated freely (360 degrees), adjustable in increments of 0.1 degree, freely adjustable XY offset
Scanning Field	20 mm field diagonal (max. 18 mm for Axio Examiner) in the intermediate image plane, with full pupil illumination
Pinholes	Master pinhole with preset size and position; can be adjusted as desired for multitasking and short wavelengths (such as 405 nm)
Beam Path	Exchangeable Twin Gate beamsplitter with up to 100 combinations of excitation wavelengths and outstanding laser line suppression; manual interface port for two external detection modules (such as NIR, BiG.2, Airyscan 2, third party detectors), internal detection with spectral signal separation and signal recycling loop for compensation of polarization effects

Detection Options

Detectors	1, 4 or 32 GaAsP PMT combined with 2 multialkali PMT internal spectral detection channels (QE 45 % typical for GaAsP) Additional Detection: 2ch NIR (GaAs and NIR GaAsP) detection or 2ch BiG.2 (UV-Vis GaAsP) detection Airyscan 2 detector (32 channels GaAsP), delivers resolution up to 120 nm lateral, 350 nm axial; Multiplex resolution: 140 / 160 nm lateral, 450 nm axial Up to 12 non-descanned detection channels (PMT and/or GaAsP) depending on microscope stand Transmitted light detector (PMT)
Spectral Detection	3, 6, or 34 + 2 NIR simultaneous, confocal reflected-light channels, GaAs, GaAsP (UV-Vis and NIR) and multialkali PMT based; freely adjustable spectral detection area (resolution down to 3 nm)
Data Depth	8 bit or 16 bit available; up to 34 channels simultaneously detectable
Real-time Electronics	Microscope, laser, scanning module and additional accessory control; data acquisition and synchronization management through real-time electronics; oversampling read-out logic; ability to evaluate data online during image acquisition

Technical Specifications

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › **Technology and Details**
- › Service

ZEN Microscopy Software

System Configurations	Workspace to conveniently configure all of the motorized functions of the scanning module, laser and microscope; Save and restore application configurations as experiment settings or use acquired images (Reuse)
Maintenance and calibration Tools	Software tools and wizards to automatically test and calibrate the system
Recording Modes, Smart Setup	Spot, Line/Spline, Frame, Tiles, Z Stack, Lambda Stack, Time Series and all combinations (XYZ, lambda, t), online calculation and visualization of, average and summation (by line/image, adjustable), Step Scan (for higher image frame rates); Quick set up of imaging conditions using Smart Setup by simply selecting the labelling dye
Crop Function	Easily select the scanning area by defining simultaneously zoom, offset and rotation
Real ROI Scan, Line and Spline curve Scan	Scan multiple ROIs (regions of interest) as desired and pixel-by-pixel laser blanking; Scan along a freely defined line
ROI Bleaching	Localized bleaching in multiple bleach ROIs for applications such as FRAP (fluorescence recovery after photobleaching) or uncaging; Use a speed or z-position different from imaging settings, use of different laser lines for different ROIs
Multitracking	Rapidly change excitation lines when recording multiple fluorescences for the purpose of minimizing signal crosstalk and increasing dynamic range
Multiplex Mode	Multiplex mode scan with 4x or 8x parallelisation in Y-direction, detection by Airyscan 2
Lambda Scan	Parallel or sequential acquisition of multidimensional images with spectral information for every pixel
Linear Unmixing	Acquisition of crosstalk-free, multiple fluorescence images using simultaneous excitation; Online or offline and automatic or interactive unmixing; Advanced unmixing logic with indication of reliability
Visualization	2D (XY); Split (XY-ch); Gallery (XY-ch, XY-Z), Orthogonal (XY, XZ, YZ) with adjustable cut lines, maximum intensity projection and 3D distance measurement; 2.5D viewing with various rendering options and animations; Histogram settings using channel specific brightness, gamma and contrast; color table selection and modification (LUT), various annotations
Image Analysis and Operations	Colocalization and histogram analysis with individual parameters, number & brightness analysis, profile measurement along user-defined lines, measurement of lengths, angles, areas, intensities and much more; operations: addition, subtraction, multiplication, division, ratio, shift, filters (low-pass, median, high-pass, etc., also user-definable)
Image Management	Features for managing images and the corresponding imaging parameters
3Dxl Viewer powered by Arivis	Rapid 3D and 4D reconstructions and animations (available modes: transparency, volume, maximum intensity projection, surface and mixed rendering with clipping planes)

Technical Specifications

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › **Technology and Details**
- › Service

Optional Software	
Direct Processing	Processing of large data during acquisition by streaming, including e.g. Airyscan, Denoising and Deblurring; analysis and storage on second PC
Deconvolution	3D, GPU based Cuda image restoration based on calculated point-spread functions (modes: nearest neighbor, maximum likelihood, constrained iterative)
HDR	Imaging mode: High Dynamic Range, improvement of the dynamic signal range by combination of multiple images with ramped signal
Physiology (Dynamics)	Comprehensive evaluation software for online and offline ratio imaging with various pre-defined formulas
FRET	Acquisition of FRET (Förster resonance energy transfer) image data with subsequent evaluation; Acceptor Photobleaching and Sensitized Emission methods supported
FRAP Efficiency Analysis	Acquisition of FRAP (fluorescence recovery after photobleaching) experiments with subsequent evaluation of intensity kinetics
RICS Image Correlation	Single molecule imaging and analysis using multialkali or GaAsP PMT detectors (publ. v. Gratton)
Experiment Designer	Definition of customized imaging configurations and procedures
Open Application Development (Macro Environment)	Python scripting interface for automation & customization; experimental feedback for smart experiments and open interface to third party software (e.g. ImageJ)
ZEN Connect and ZEN Connect 2D/3D Add-ons	Exchange and alignment of image data from multiple image acquisition systems in 2D and 3D enabling correlative workflows
ZEN Intellesis	Image analysis and structure detection via computational self learning technology
FCS/FCCS	Fluorescence Correlation and Cross Correlation Spectroscopy for analysis of single molecule dynamics, concentration and number
AI Sample Finder, Sample Navigator (requires additional HW)	Easy to perform sample overview scan with autofocus function using AxioCam or transmitted fluorescence with T-PMT (Finder requires Axio Observer)
Guided Acquisition	Automated and targeted acquisition of objects of interest
Tiles & Positions	Scanning of predefined sample areas (tiles) and/or position lists
3Dxl Plus	Combine 2D and 3D visualization in one screen
BioApps	Easy-to-use and modular image analysis for common assays
Airyscan RAW data	Optional export of complete Airyscan single channel data and the sheppard sum for external processing, e.g. correlations, deconvolution, AI etc.

Technical Specifications

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › **Technology and Details**
- › Service

Lasers

Laser RGB (445, 488, 514, 543, 561, 594, 639 nm)

Single-mode polarization preserving fiber

Laser beam attenuation for all lasers by VIS-AOTF

Diode Laser 445 nm (30 mW nominal power; 7.5 mW ex fiber)

Diode Laser 488 nm (30 mW nominal power; 10 mW ex fiber)

Diode Laser 514 nm (30 mW nominal power; 10 mW ex fiber)

DPSS Laser 543 nm (25 mW nominal power; 10 mW ex fiber)

DPSS Laser 561 nm (25 mW nominal power; 10 mW ex fiber)

DPSS Laser 594 nm (8 mW nominal power; 2.5 mW ex fiber)

Diode Laser 639 nm (25 mW nominal power; 7.5 mW ex fiber)

Laser V and NIR (405 and 730 nm), direct modulated

Single-mode polarization preserving fiber

Diode Laser 405 nm (30 mW nominal power; 14 mW ex fiber)

Diode Laser 730 nm (10 mW nominal power; 9.5 mW ex fiber)

Power Requirements

LSM 980 has a main power supply cord and plug, either NEMA L5-15 (100V – 125V) 2pol (15A) + PE or CEE blue (200 – 230V) 2pol (16A) + PE.

Line Voltage 1/N/PE 230 V AC (±10 %)

1/N/PE 120 V AC (±10 %)

Line Frequency 50...60 Hz

50...60 Hz

ZEISS LSM 980 incl. VIS Laser

Max. Current 7 A at 230 V

13 A at 120 V

Heat emission without Multiphoton Laser 1,500 W max.

1,500 W max.

Power Consumption 1,600 VA max.

1,600 VA max.

Multiphoton Laser

Power consumption and heat emission varies depending on type of laser. See data sheet of laser from laser supplier for further information.

EMC test

according to DIN EN 61326-1

– Emitted interference according to CISPR 11 / DIN EN 55011

– Interference immunity as specified in Table 2 (industrial)

Nominal power = power level of laser itself, leaving out necessary tolerances and losses due to laser steering and stability requirements



Technical Specifications

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › **Technology and Details**
- › Service

Environmental Requirements

For operation the system has to be placed in a closed room.

1. Operation, specified performance	T = 22 °C ± 3 °C without interruption (24 h a day independently whether system is operated or switched-off) It has to be ensured that the air-flow of the air-conditioning is not directed at the system.
2. Operation, reduced performance	T = 15 °C to 35 °C, any conditions different from item 1. and 5.
3. Storage, less than 16 h	T = -20 °C to 55 °C
4. Temperature gradient	± 0.5 °C/h
5. Warm up time	1 h, for high-precision and/or long-term measurements ≥ 3 h
6. Temperature gradient and range for continuous long-term image acquisition	± 1.5 °C/12 h
7. Relative humidity	< 65 %
8. Operation altitude	max. 2,000 m
9. Loss of heat (without Multiphoton Laser)	1.5 kW
10. Vibrations under operation conditions (with system table)	Vibration Class 12.5 µm/s VC-C (IEST RP 12 and ISO 10811)
11. Shipping shock (LSM 980 box)	< 10 g

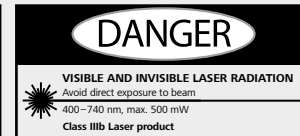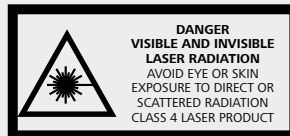

LSM 980 meets the requirements according to IEC 60825-1:2014

- › In Brief
- › The Advantages
- › The Applications
- › The System
- › Technology and Details
- › **Service**

ZEISS Service – Your Partner at All Times

Your microscope system from ZEISS is one of your most important tools. For over 170 years, the ZEISS brand and our experience have stood for reliable equipment with a long life in the field of microscopy. You can count on superior service and support - before and after installation. Our skilled ZEISS service team makes sure that your microscope is always ready for use.

Procurement

- Lab Planning & Construction Site Management
- Site Inspection & Environmental Analysis
- GMP-Qualification IQ/OQ
- Installation & Handover
- IT Integration Support
- Startup Training

Operation


- Predictive Service Remote Monitoring
- Inspection & Preventive Maintenance
 - Software Maintenance Agreements
 - Operation & Application Training
 - Expert Phone & Remote Support
- Protect Service Agreements
 - Metrological Calibration
 - Instrument Relocation
 - Consumables
 - Repairs

New Investment

- Decommissioning
- Trade In

Retrofit

- Customized Engineering
- Upgrades & Modernization
- Customized Workflows via APEER

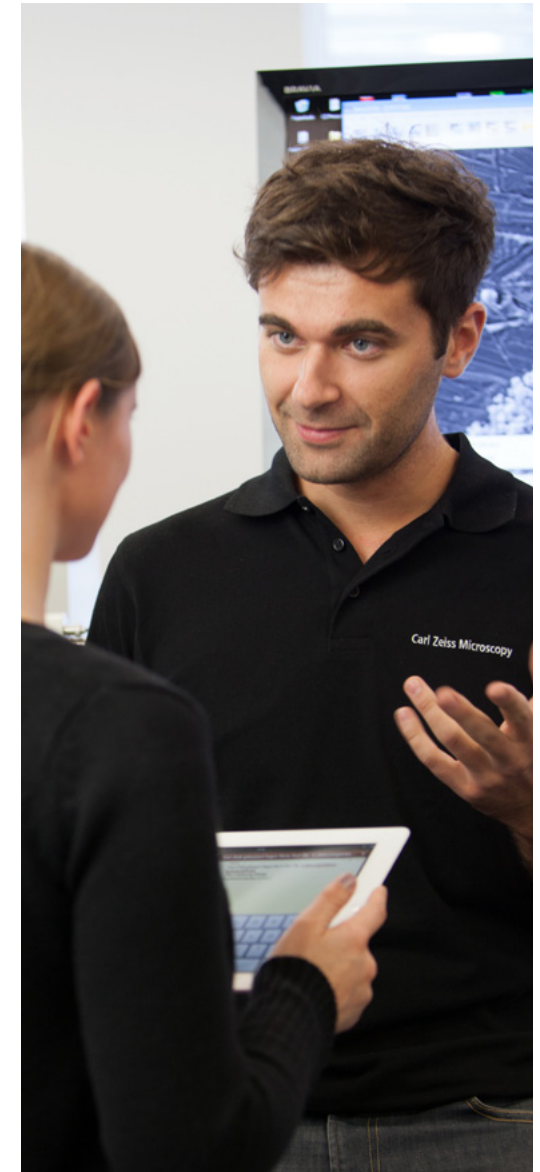

Please note: Availability of services depends on product line and location

>> www.zeiss.com/microservice



Carl Zeiss Microscopy GmbH

07745 Jena, Germany
microscopy@zeiss.com
www.zeiss.com/lsm980

